MODA sofa

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User Manual
OS 3008 使用手册

MODA sofa

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User Manual
OS 3008 使用手册

CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Controller Diagram	07
Button Name and Function	08
Product Installation	09
Product Operation Guide	12
Display	13
Product Maintenance	14
Trouble Shooting	15
Product Specification	15
Contact us	16

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Place the appliance in a safe location.
- Do not expose the appliance to a heater or direct sunlight.
- This product is only intended for personal use only. Commercial use is prohibited.

ii. Health precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease.
- Pregnancy.
- Any health concerns.

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Before sleeping.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

iv. Caution and warning

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by OGAWA.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.
- Do not insert your hands or any foreign matters into the massager's vents and gaps.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- This product is only intended for personal use only. Commercial use is prohibited.

Safety Instructions

- Do not cover the appliance during operation. As this may result in malfunction or fire.
- Do not crush the appliance and avoid any sharp folds.


V. Note

- This appliance can be used by children at the age of 8 or above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under guidance or if they have been informed of the risks involved and instructions concerning use of the appliance in a safe way.
- Children shall not play on the appliance.
- Cleaning and maintenance shall not be performed by children without supervision.
- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- Extensive usage could lead to excessive heating and a shorter product lifespan. Should this occur, discontinue the use and allow the unit to cool before operating again.
- It is normal to hear the noises during use, because they are resulted from vibration of product.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.
- If damaged, the supply cord must be repaired by the manufacturer, after-sales service provider or qualified persons in similar field so as to avoid a hazard.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Maximum user weight is 120Kg.


Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

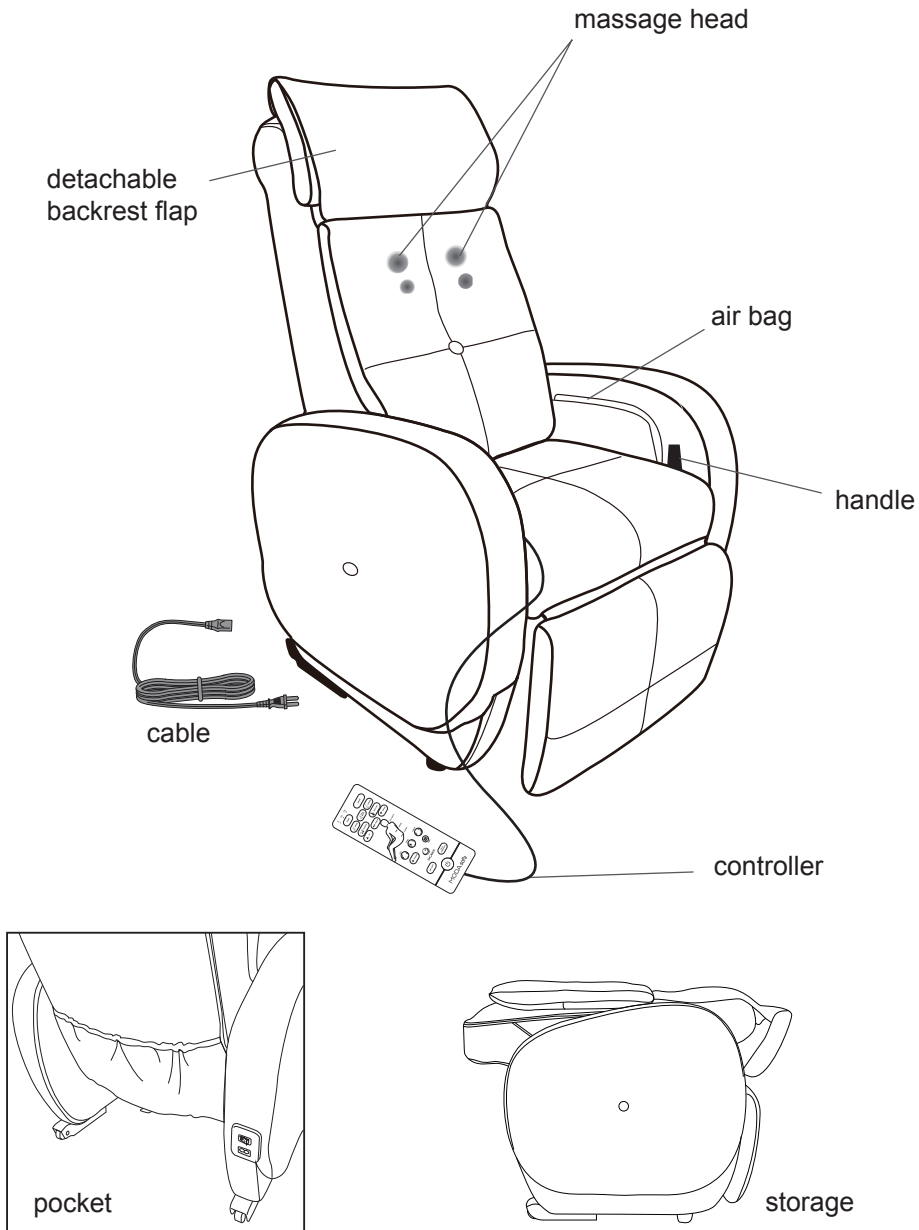Controller Diagram

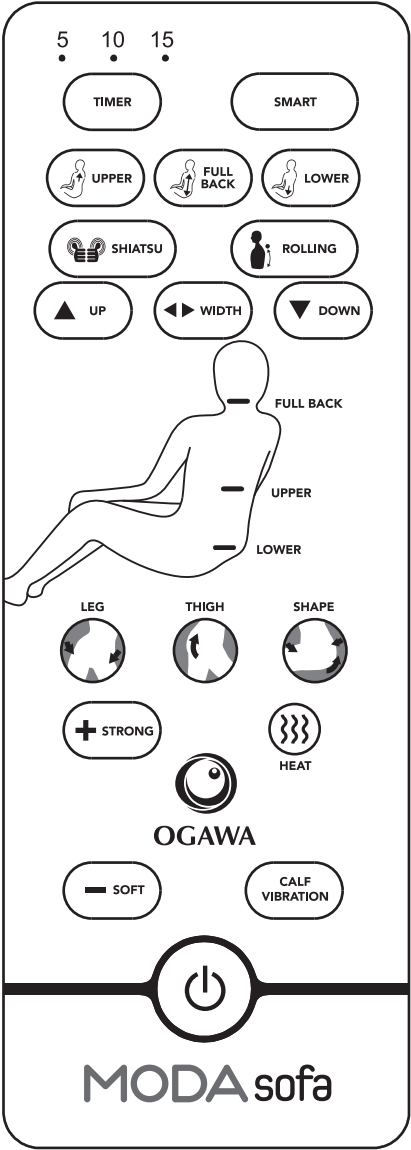
Button Name and Function

1. Timer - use to select massage time
2. Smart - for back and thigh massage
3. Area - use to select massage area
4. Shiatsu - use for back shiatsu massage
5. Rolling - for back and thigh massage
6. Width - use to select massage head width
7. Up/Down - use to select massage position
8. Let - Leg massage program
9. Thigh - thigh massage program
10. Shape - leg and thigh massage program
11. Intensity - air bag pressure selection
12. Heat - controller button for heat at back
13. Vibration - controller button for vibration
14. Power - button to ON/OFF machine


Product Installation

Step 1. Lift up the backrest, Make sure the backrest was lock
(locking sound "ka")

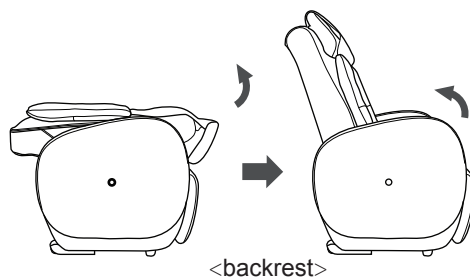
Step 2. Assemble the connection by follow
the indication

Caution: Unproper connection will cause
machine unoperated.

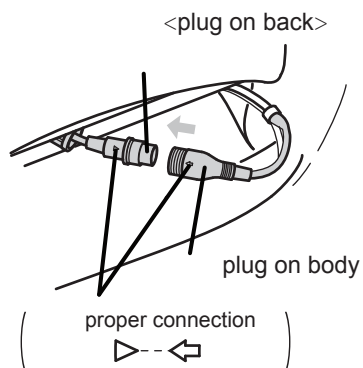
Step 3. Plug in socket.
as picture

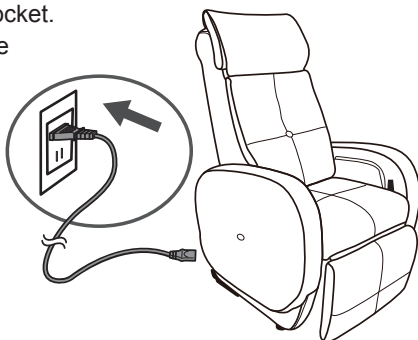Back rest storage

Caution: Make sure the massage head is return to bottom
positon before dismantel. Remove the plug.

1. Remove the connection
2. Press the button to release
the backrest as picture

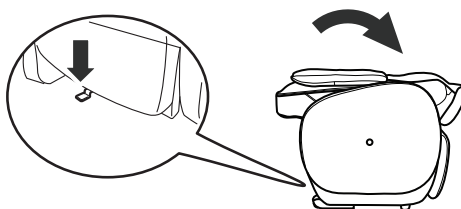Product Operation Guide

1. [Power]

Insert power plug into electric outlet. Turn on the power. If the roller are not at bottom position, the light of the [Power] button will flash and the rollers will return to the bottom position. Press [Power] button, light will on and start countdown of 15 minute. Press [Power] second time, the unit will stop functioning and the light is off.

2. [Timer]

Each massage is set for 15 minute, each time button was press, the time will change in cycle 5-10-15. The light in each location will on.

3. [SMART]

Press [SAMRT] button, whole body, shiatsu, heat , thigh and leg massage. air pressure in low level for 15 minute. [Heat] and [Intensity] can be press to switch off or change level. If other button is press(except width), [SMART] will be off and proceed to the related function. Press [SMART] again to turn off this program.

Back massage

1 . [Shiatsu]

Press [Shiatsu] to switch on the shiatsu massage. Press again the stop shiatsu massage.

2. [Up / Down]

Long pres [Up] button, massage head will move upward and stop when button is release.

Long pres [Down] button, massage head will move upward and stop when button is release.

3. [Rolling]

Press [Rolling] button, rolling massage will activate. Press [Rolling] button again, rolling massage will stop.

4. [Width]

Width only will be function under rolling mode.

During rolling, each time when [Width] is press, massage head turn to provide different width.

Product Operation Guide

5. [Full Back] , [Upper] , [Lower]

[Full Back],[Upper],[Lower] button is to select the massage area. If none of the button is press, machine default as full back, if none of [Shiatsu], [Rolling] button is press, machine will default as shiatsu. Press [Full Back] button, button light will on and massage head move all location. Press [Upper] button, button light will on and massage head move only upper half. Press [Lower] button, button light will on and massage head move only lower half.

6. [Heating]

Heating only function when shiatsu, rolling is activated. Press [Heating] button, heating will start. Press again to stop it.

7. [Vibration]

Press [Vibration] button, vibration will activated. Press again it will change to wave type vibration. Press again will change to motion type vibration; Press again to stop it.

Air pressure massage

1. [Leg]

Press [Leg] button, button light will on, masage on leg and hip.

2. [Thigh]

Press [Thigh] button, button light will on, masage on thigh.

3. [Shape]

Press [Shape] button, button light will on, air bag on hip and thigh will exchange massage.

4. [+Strong] [-Soft]

Use to select air pressure intensity. Default as low pressure. Last setting will be memory if [Power] button not press to off.


Caution:

- Please use backrest flap if found the massage is too strong.
- For new user, please select low air intensity.
- After use, press power button to off the machine. Once the massage head had back to the origin.

Product Operation Guide

5. Backrest angle control

Backrest incline decline angle can be adjust.

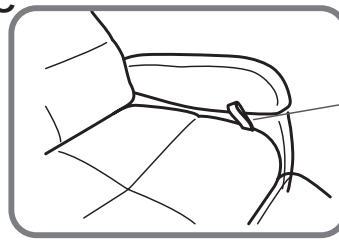handle

A. Decline:


Before decline, please make sure no disturbance behind chair. Pull up handle, lay down backrest.

B. Incline:

Pull up handle, backrest will back to the origin position.

6. Relocation:

Design with roller at rear bottom. To shift, press down chair as picture and pull.

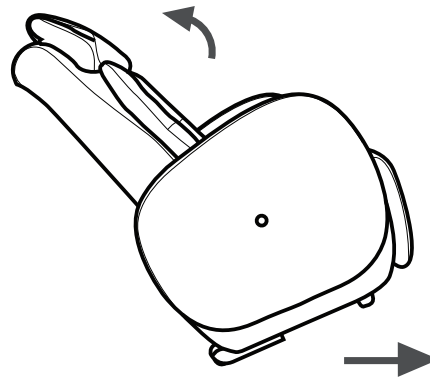Caution:

- Please use backrest flap if found the massage is too strong.
 - For new user, please select low air intensity.
 - After use, press power button to off the machine.
- Once the massage head has back to the origin, remove the power plug from machine.

Display

1. Shiatsu and rolling massaging programs can relax your back.
2. Fixed-position massaging function provides Shiatsu of specific position.
3. Leg cushion vibration relieves the leg tension.
4. Heat function can increase the Shiatsu and Rolling effect and improve blood circulation.
5. Easy to use remote controller.
6. Auto shut-off after 15 minutes.
7. LED light indicators.

Effects by using the appliance

1. Relieves the muscle ache and pain.
2. Relieves the tension and reduce the stress.
3. Improves blood circulation.

When you do not use the product, press the [Power] button, and please don't cut off the power immediately until the rollers have resumed to the bottom and the indicator light is off.

Press the [Power] button after use and the massager will stop operating. Massage roller will resume to their original position. Unplug the appliance from power outlet.

Product Maintenance

TO STORE

- Keep massager in a safe, dry and cool place.
- Keep away the massager from sharp edges or pointy objects that might cut or puncture the product's fabric surface.
- Do not wrap the power cord around the unit.

TO CLEAN

- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Keep the appliance away from water or any liquid.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish or paint thinner to clean the product.

CAUTION

- Do not lean on the dual massage mechanism in the cushion when adjusting your body position to avoid tears on the products cushion surface.
- Do not jam or force any part of your body against the moving massage mechanism.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.
- Do not attempt to repair the product, please contact 'OGAWA Service Center' for any servicing required.

Troubleshooting

Symptoms	Solutions
Not operating.	<ul style="list-style-type: none"> • Use the product when electricity is restored. • Reconnect the power cord. • Turn on the power switch.
No heating after pressing the heat switch.	<ul style="list-style-type: none"> • Press the 'shiatsu' switch before turning the 'Heat' switch. • Press the 'Rolling' switch before turning on the 'Heat' switch.
During operation, the controller is working but the massager is not working.	<ul style="list-style-type: none"> • The massager could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and stop using the massager for 30 to 50 minutes to facilitate heat dissipation of the engine.
Back massage function does not work.	<ul style="list-style-type: none"> • Turn on the power switch. • Reconnect the power cord.
Leg vibration does not work.	<ul style="list-style-type: none"> • Turn on the 'Vibration' button. • Reconnect the power cord.
Sudden stop.	<ul style="list-style-type: none"> • Use the product when electricity is restored. • Reconnect the power cord.

Product Specifications

Product Name	MODA sofa
Model number of the product	OS 3008
Power source	220-240V
Frequency	50~60Hz
Power consumption	72W
Dimension (Length/width/height)	Decline:700*720*110 (mm)
	Incline: 1300*720*950 (mm)
Weight	32.6KG

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F, YHC Tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel : +852 3465 6262
Fax: +852 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM
OGAWA Vietnam
171 Nguyen Thi Thap Street
Tan Phu Ward
Dist 7
Ho Chi Minh City, Vietnam
Tel : +84 8 3775 3093 / 3775 3094
Fax: +84 8 3775 3095

INDONESIA
PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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