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
210mm

IntelLAB
bluetooth body scale

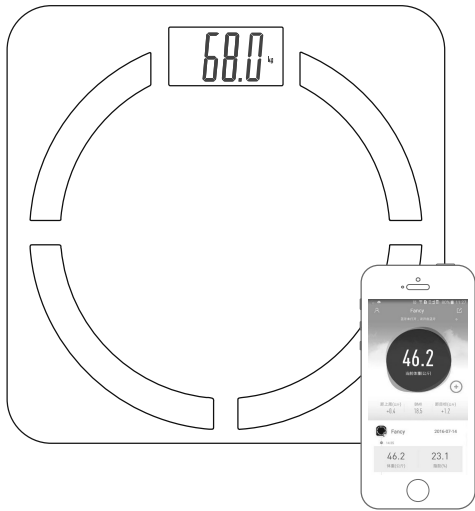User Manual
AC 1202 使用手册

Contents

Instruction

- Product instruction 1
- Safety warning 1
- Main parts 1
- Battery installation 2

Measuring

- APP Instruction 3
- APP Installation 4
- Registration 5
- measuring 6
- Pairing / unbinding 6-7
- Measure data and report 7-8
- Family members adding and recognition 8

Features

- recycling 10
- Maintain 10
- Warranty Claim 11
- Contact Us 12

Instruction


Product instruction:

Bluetooth scale helps you to store and keep track of your weight on mobile devices by OKOK app. This scale can transmit data of your Body fat, body hydration, body muscle, body bone, metabolism, BMI, Visceral Fat, Protein, Obesity Degree, Body Age, Weight without fat etc.to your smart devices via bluetooth.

Safety warning :


1. Those who use heart pacemaker or with other implanted medical device should not use the device. That's because Bio-electric current through the body when measuring. Then the medical devices may not work normally.
2. The Fat%,hydration% content evaluated by the scale is only for reference(not for medical purpose). If your Fat%,hydration% content is over/below normal level ,please consult your doctor for more advice.
3. People with below conditions are not advisable using the device: children under 10, elder above 80, and pregnant women.

Main parts:

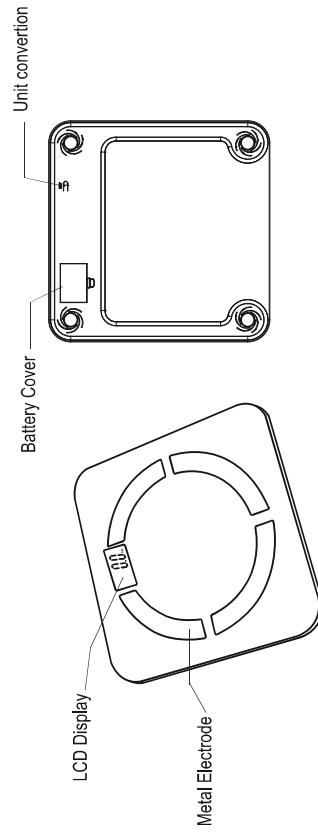Battery installation

1. Open the battery cover on the backside of the product.
2. Install the batteries and ensure + /- terminals are correct.
3. Close the battery door.

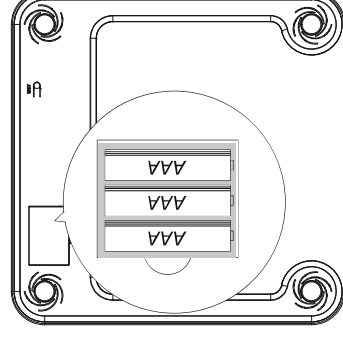
Under below condition, please change new batteries

- Display 
- Dimmer display
- Nothing display




※ Note:

- Please use corresponding battery.
- Please remove battery if it is not being used for a long time.

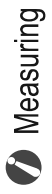
APP Instruction


Main page overview

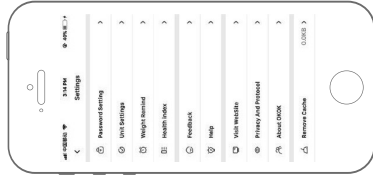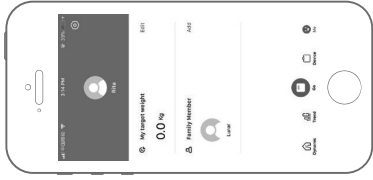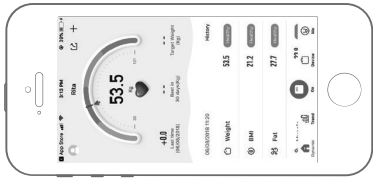8 users' memory

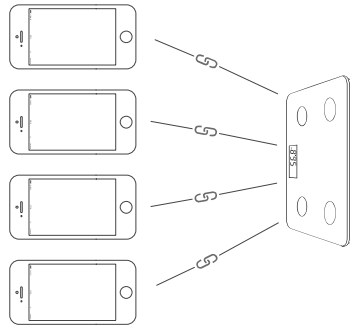
Detailed measure data

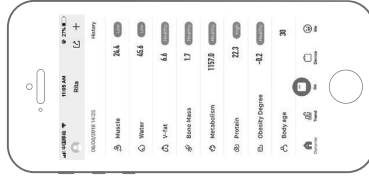
- Body fat
- Body hydration
- Body muscle
- Body bone
- Metabolism
- 8 users
- BMI
- Visceral Fat
- Protein
- Obesity Degree
- Body Age
- Weight without Fat

APP Installation

1. Compatible with the following mobile OSes:

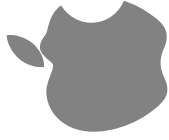IOS 7.0

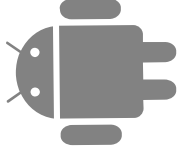Android 4.3

Bluetooth 4.0

2. Turn on your mobile phone's bluetooth and WIFI or 3G/4G.

3. Scan the QR below to install OKOK app.

also you can download OKOK app from the App store(for IOS only).


Registration


Login


profile setting

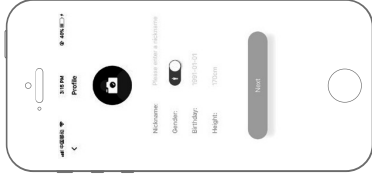
Target weight setting

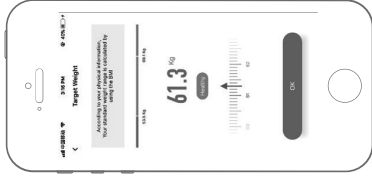measuring

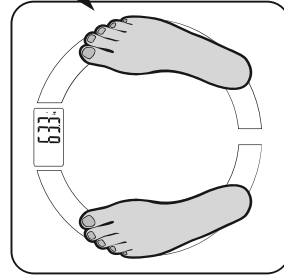Make sure good contact with the metal electrodes


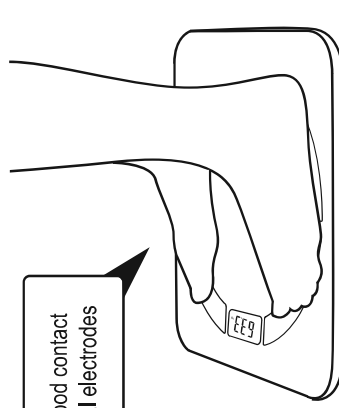
Place the scale on a hard, flat platform

Tips:

1. Step on the scale with bear feet.


Pairing

Click "Go" on the main page to finish the binding.

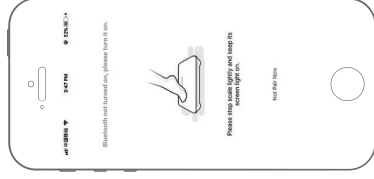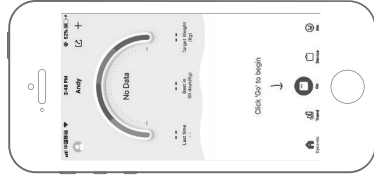unbinding

Press "Device" on the main page to finish unbinding

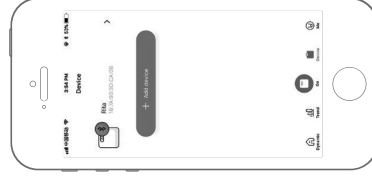Measure data and report

Showing weight

Showing all measure values

Health data analyze & health tips


In the health score page, tap on every measurement, the analysis and health tip will come out

Body fat, body hydration, body muscle, body bone, metabolism, BMI, Visceral Fat, Protein, Obesity Degree, Body Age, Weight without fat


Family members adding and recognition

Press "Me" on main page, to add more family members

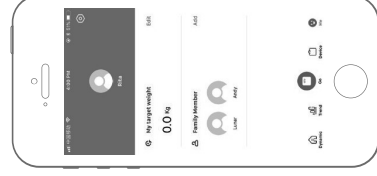A total of 8 family members can be added.
Using a 3-kg difference in weight, the scale automatically recognise the family member.

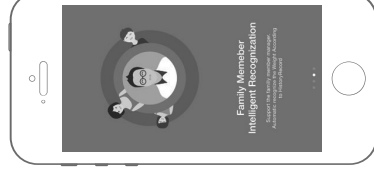
Features

Item	Parameter
Capacity	180kg/397lb/28st:5lb
Division	0.1KG/0.2lb/1lb
Unit	kg; lb; st:lb
High fat indication	ErrH
Low fat indication	ErrL
Height range	100 ~ 240/cm;3' 3 ~ 7' 10/ft inch
Measuring body fat range	5% ~ 50%
Measuring water range	10% ~ 80%
Measuring bone range	0.5kg ~ 8.0kg
Measuring muscle range	10% ~ 80%
Auto identification range	3kg


recycling

When no longer functional an electronic product should be disposed of with the least possible environmental impact. The device should be disposed of according to your local municipal regulations, but in most cases you can take the product to your local recycling centre.

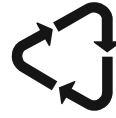
Product maintain


Please follow below instruction to maintain the product, so that to keep it in good condition.

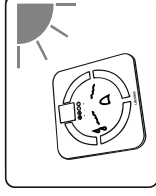Avoid to use the scale in humid, overheat or overcold environment. (temperature range from $-20\sim 60^{\circ}\text{C}$)

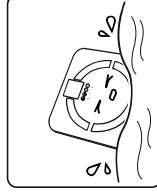Do not put the product in water or wash it with water, to clean the scale with a slightly damp cloth when it is dirty.

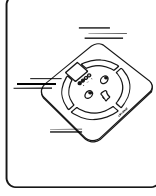Do not drop, shock or strike the scale as it is a precise instrument.

Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the message unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Contact Us

MALAYSIA

Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA

OGAWA Health Care International (HK) Ltd
31/F, YHC Tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel : +852 3465 6262
Fax: +852 2410 0119

SINGAPORE

OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM

OGAWA Vietnam Co Ltd
2nd floor, Petroland Tower,
12 Tan Trao, Tan Phu ward
District 7
Ho Chi Minh
Tel : +84 8 541 33 222 / 54133 223
Fax: +84 8 5416 1736

INDONESIA

PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES

O-Healthcare Solution Phil.Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / + 632 844 6851
Fax: +632 812 0421

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O-Healthcare Solution Phil. Inc.
Thailand . Korea . Australia . Dubai . K.S.A . Nigeria . Egypt . Lebanon . Denmark
Sweden . U.S.A