
ACU THERAPY

Acupressure + Air Compression
Foot Massager with Heat

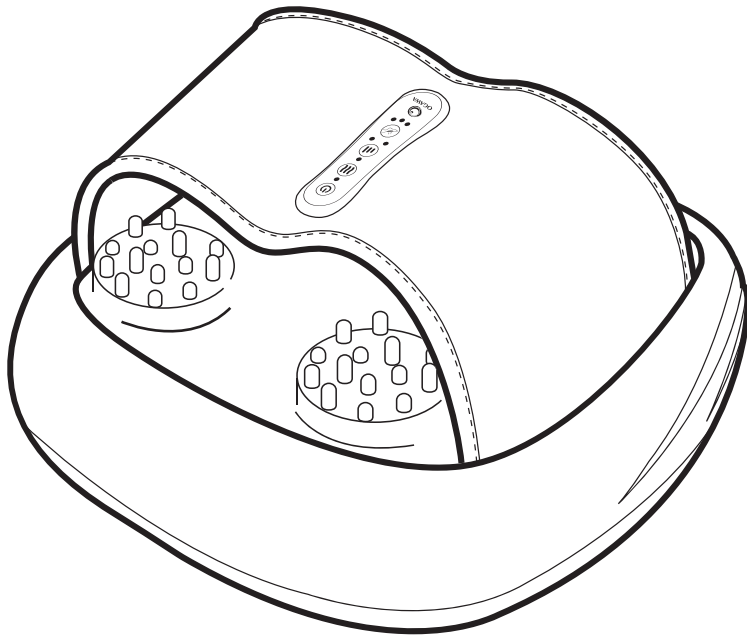User Manual
OF 0609 使用手册

ACU THERAPY

Acupressure + Air Compression Foot Massager with Heat

User Manual

All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Controller Diagram	06
Product Operation Guide	07
Product Maintenance	08
Troubleshooting	09
Product Specifications	10
Contact Us	11

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not use the appliance on the bed.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink. Place the appliance in a safe location.
- Do not sit on the appliance.
- Do not expose the appliance to a heater or direct sunlight.

ii. Health Precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease
- Pregnancy
- Any health concerns

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Do not use this product before sleep.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

iv. Caution and warning

- Unplug the appliance from the outlet when not in use
- Do not use any attachments that are not recommended by OGAWA
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.

Safety Instructions

- Do not drop or insert any object into any opening.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- Do not lift or carry this appliance by the power cord or any cord.
- This appliance is designed only for personal use, and not for professional purposes.
- Do not cover the appliance during operation.
- Do not crush the appliance and avoid any sharp folds.



v. Note

- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- The appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. They need to be supervised and instructed about its usage by a person responsible for their safety.
- Extensive usage could lead to excessive heating and a shorter product lifespan. Should this occur, discontinue the use and allow the unit to cool before operating again.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.

Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

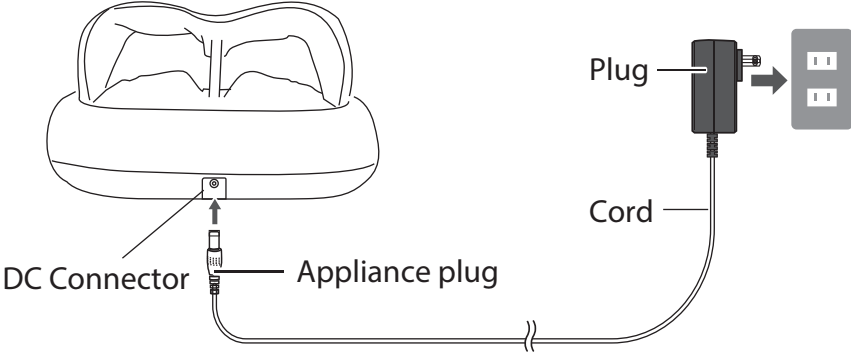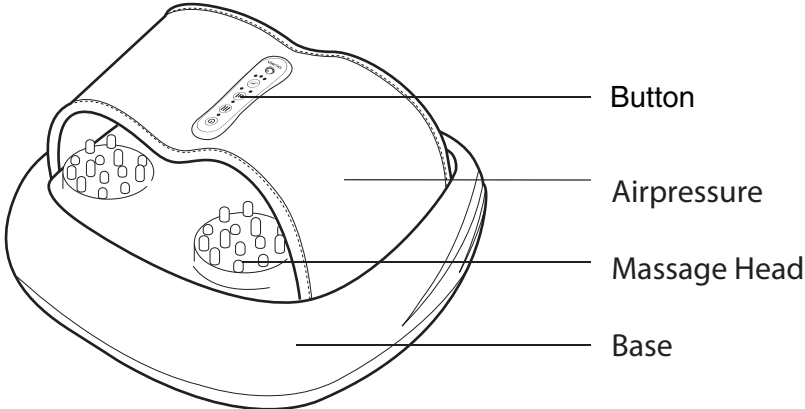Product Operation Guide


i. Preparation before use

1. Before using the appliance
 - Check if there is any damage to the appliance.
 - Check the appliance is placed in the proper position.

ii. Operation guide

1. **(Power ON/OFF):** Press this button, machine power-on. All keys are effective. Automatically setting time is 15 mins. The Acupuncture and Air pressure function turn on when power on, and the air-pressure of the feet is opened by default in weak intensity, Acupuncture speed defaults to slow speed. Press this button again, the whole machine is shut down, and all massages are stopped and the system is in standby mode. When time's up, machine stop massaging, enter standby state.
2. **Speed:** It is defaulted to be in Slow speed when power on, press the button to adjust the speed as you need. Speed switches every time when press button, shows by 2 LED indicator lights, respectively indicate the Slow and Fast speed.
3. **Heating:** It is defaulted to open heating function in startup state. Press button to close Led light and heating function. Followed by recycling.
4. **Air pressure:** It is defaulted to be in Weak Intensity on power-on state, Press button to switch intensity level to be in Medium, corresponding indicator light on.

1. Power button ON/OFF
2. Heating button: ON/OFF
3. Massage speed (Low / High)
4. Air Compression (Low / Medium / High)

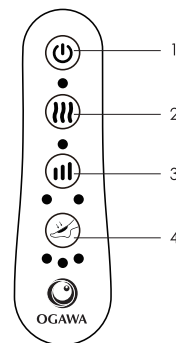Product Maintenance

TO STORE

- Place the massager in this box or in a safe, dry and cool place. Avoid any contact with sharp edges or pointed objects that might cut or puncture the surface.
- Please wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

TO CLEAN

- Close the power button and unplug the unit before cleaning.
- Clean only with a soft, slightly damp sponge.
- The attachable cloth can be used to clean the appliance.

CAUTION

- Avoid using any corrosive detergent, such as Benzene or diluents to clean the appliance.

MAINTENANCE

- Keep the appliance clean and store it away from wet conditions.
- If the appliance is not used for a long time, please use a dust cover (please buy it by yourself) to cover the appliance, to avoid any dust on the product.

Operations for Use

Symptom	Cause	Solution
The product does not work	The plug is not inserted properly	Insert the plug properly and safely
	The power key is not started	Start the power key
	The product will be turned off automatically after it finishes a 15 minutes' massage	If you want massage for another time, press power key again
The product is hot or cannot be used suddenly in midway	You have performed several groups of 15 minutes' massage already.	Take off the power plug and let it cool. Always take off the power plug after using it.

If you still have any problem after using the methods above, please contact your service center .

Product Specifications

Product Name	Acu Therapy
Model number of the product	OF 0609
Power source	DC 12V --- 2.5A
Power consumption	30W
Timer	15 minutes
Dimensions (Length/width/height)	L330 x W160 x H360 (mm)
Weight	3.0 ± 0.1kg

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F,YHC Tower,1Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel:+852 3465 6262
Fax:+852 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampine Industrial Avenue 5,
T5 @ Tampines
Singapore 528623

Tel:+065-6846 8588
Fax:+065-6846 7677

VIETNAM
OGAWA Vietnam Co Ltd
2nd floor,Petroland Tower,
12 Tan Trao,Tan Phu ward
District 7
Ho Chi Minh
Tel:+84 8 3775 3093/3775 3094
Fax:+84 8 3775 3095

INDONESIA
PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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