Tapping Foottee

foot reflexology *therapy*

User Manual
OF 1728 使用手册

OGAWA

Tapping Foottee

foot reflexology *therapy*

User Manual

All information is correct at the time of print

CONTENTS

| | |
|-------------------------|----|
| Safety Instructions | 03 |
| Warranty Claims | 05 |
| Main Unit & Accessories | 06 |
| Controller Diagram | 06 |
| Product Operation Guide | 07 |
| Product Maintenance | 08 |
| Troubleshooting | 09 |
| Product Specifications | 10 |
| Contact Us | 11 |

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not use the appliance on the bed.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink. Place the appliance in a safe location.
- Do not sit on the appliance.
- Do not expose the appliance to a heater or direct sunlight.

ii. Health Precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease
- Pregnancy
- Any health concerns

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Do not use this product before sleep.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

iv. Caution and warning

- Unplug the appliance from the outlet when not in use
- Do not use any attachments that are not recommended by OGAWA
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.

Safety Instructions

- Do not drop or insert any object into any opening.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- Do not lift or carry this appliance by the power cord or any cord.
- This appliance is designed only for personal use, and not for professional purposes.
- Do not cover the appliance during operation.
- Do not crush the appliance and avoid any sharp folds.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

v. Note

- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- The appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. They need to be supervised and instructed about its usage by a person responsible for their safety.
- Extensive usage could lead to excessive heating and a shorter product lifespan. Should this occur, discontinue the use and allow the unit to cool before operating again.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.


Safety Instructions

- The appliance is only to be used with the power supply provided.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.


Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

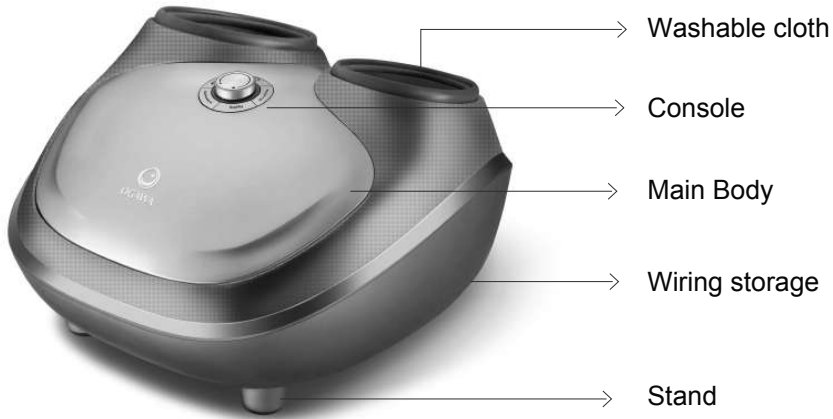Controller Diagram

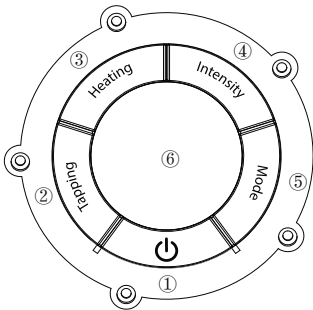- 1) POWER Button-Control supply of the appliance, When button is press, the LED light will on.
- 2) Tapping Button-3 selection mode
Mode 1: green LED is on
Mode 2: blue LED is on
Mode 3: red LED is on
- 3) Heating Button-Press this button to activate heat function, the red LED will on. Press again to stop the function.

4) Intensity Button-Control the massage intensity of the appliance via 3 modes:
The green LED is weak, blue LED is middle, red LED is strong.

5) The Mode Function-
Mode 1: blue LED is on
Mode 2: green LED is on
Mode 3: red LED is on

6) Tapping speed controller- CW is to increase the speed. CCW is to decrease the speed.

Product Operation Guide

i. Preparation before use

1. Before using the appliance
 - Check if there is any damage to the appliance.
 - Check the appliance is placed in the proper position.

ii. Operation guide

1. Untie the cord and plug in the appliance into 220-240V 50/60Hz socket. the buzzer will ring, the appliance will enter a standby state.
2. Operating method
 - The appliance are in stand by mode, touch Power button, buzzer will sound, decorative light will on, product enter to on mode, at the same time appliance enter to 15 minute countdown.
 - Touch Mode button to activate Mode 1, Mode 2, Mode 3-close(at the same time shut off tapping & heating) in a cycle. Each change, the buzzer will ring and corresponding LED will turn on.
 - In auto mode, intensity is preset in week, green LED, tapping is in standby mode. white LED.
Mode 1: blue LED is on
Mode 2: green LED is on
Mode 3: red LED is on

Tapping button

- 1) Tapping function can be activate when Auto mode is off. The strength can be adjust from tapping speed controller. Touch tapping button to activate tapping function.
Mode 1: green LED is on
Mode 2: blue LED is on
Mode 3: red LED is on
- 2) In auto mode,tapping button is white LED, tapping function can be on/off independently, The strength can be adjust from tapping speed controller.
 - Heating function can be on/off independently by touching the heating button. When acivate, LED light will on and LED will off when function is off. The buzzer will sound and corresponding LED will on.
 - Intensity button--Control the massage intensity of the appliance via 3 modes.
 - The green LED is weak, blue LED is middle, red LED is strong.

Operating Guide

- Tapping speed controller
Under activating tapping function, speed controller can be use
 - 1) Turn button in CCW, tapping strength will decrease.
 - 2) Turn button in CW, tapping strength will increase.
- When the appliance is on, it will enter into 15 minute countdown. Once reach 15 minute, rall process will stop. Appliance turn into standby mode.
- When the appliance is on standby mode, unplug the cord it will cut off the current.

Product Maintenance

TO STORE

- Place the massager in this box or in a safe, dry and cool place. Avoid any contact with sharp edges or pointed objects that might cut or puncture the surface.
- Please wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

TO CLEAN

- Close the power button and unplug the unit before cleaning.
- Clean only with a soft, slightly damp sponge.
- The attachable cloth can be used to clean the appliance.

CAUTION

- Avoid using any corrosive detergent, such as Benzene or diluents to clean the appliance.
- Take care of the LOGO area in order to avoid tearing the cloth via its edges.

MAINTENANCE

- Do not discard the protective film as much as possible in order to protect the metal LOGO.
- Keep the appliance clean and store it away from wet conditions.
- If the appliance is not used for a long time, please use a dust cover (please buy it by yourself) to cover the appliance, to avoid any dust on the product.

Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

| Conditions | Solutions |
|--|--|
| Not operating | <ol style="list-style-type: none">1. Check if the plug is connected properly.2. Press the "POWER" button. |
| During operation, the appliance is not operating | <ol style="list-style-type: none">1. Check the plug is connected properly.2. Press the "POWER" button again.3. After the appliance has cooled down, press the " POWER" button again. |

Product Specifications

| Product Name | Tapping Foottee |
|----------------------------------|-------------------------|
| Model number of the product | OF 1728 |
| Power source | 220-240V |
| Power frequency | 50/60Hz |
| Power consumption | 50W |
| Timer | 15 minutes |
| Dimensions (Length/width/height) | L410 x W390 x H261 (mm) |
| Weight | 9.34kg |

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F,YHC Tower,1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel :+852 3465 6262
Fax:+852 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM
OGAWA Vietnam
171 Nguyen Thi Thap Street
Tan Phu Ward
Dist 7
Ho Chi Minh City, Vietnam
Tel : +84 8 3775 3093 / 3775 3094
Fax: +84 8 3775 3095

INDONESIA
PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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