OmKNEE

total foot reflexology *therapy*

User Manual
OF 2002A 使用手册

OGAWA

OmKNEE

total foot reflexology *therapy*

User Manual

All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	06
Product Parts	07
Product Features And Operation Illustrations	08
Product Operation Guide	09
Product Maintenance	13
Troubleshooting	14
Product Specifications	15
Contact Us	16

Thank you for purchasing OGAWA OmKnee leg massager. Please read the Safety Precautions carefully before using the leg massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the degree of danger. This book should be kept available for future reference after reading.

Safety Instructions

a) Important safety instructions

- To ensure a safe and proper operation, please consult the instruction manual carefully before using this product.
- Use this product only as described in this manual.
- Always unplug this product from the electrical outlet immediately after use and before cleaning.
- Do not use this product if you have sensitive skin, neuropathy or poor circulation. Don't let children or people with reduced physical, sensory or mental disabilities use without proper supervision. Improper use may be dangerous.
- Do not operate this product under a blanket or when covered. Over-heating can occur and cause fire, electric shock or injury.
- Do not cover this product in order to avoid electrical faults or damages due to excessively high motor temperature.
- Stop the massage immediately if you feel any discomfort.
- Never operate this product if it has been dropped, has a damaged cord or plug; experiences a malfunction, or if it has been dropped in water. Please contact an authorized service centre to do the necessary check and repair.
- Do not use the components that are not recommended by OGAWA.
- Do not insert any objects into any openings in the product.
- If you want to cut off the power supply, please make sure you have turned off all the operating controls, as well as remove the power plug from the socket.
- Please do not use this product as a medical device. If you need a professional diagnosis or treatment, please consult a professional medical person.

b) Operating Environment Precautions

- Do not use this massager in humid or extremely dusty environments to avoid malfunction.
- Do not use this massager at a temperature above 40 degree centigrade.
- Do not use this massager outdoors.
- Avoid placing or operating this massager in direct sunlight or near a heating appliance, stove.

Safety Instructions

c) User precautions

- Make sure that this product is far away from any hot surfaces.
- Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, back pain, spinal scoliosis, hypertension, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc.
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.

Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

d) Precautions for power supply

- Do not turn on the switch when your hands are wet.
- Do not immerse in or spray water or liquid on the control power switch.

e) Cord/Plug Precautions

- Do not connect or remove the plug with wet hands. Otherwise you will suffer an electric shock or a product malfunction.
- Do not pull the wires when you pull out the power plug. Please hold the plug itself.
- As the default operating time is 15 minutes, all the functions will stop working after 15 minutes. Please unplug this product from the electrical outlet immediately after use.
- Keep the cord away from hot surfaces.
- Do not place the power plug under the product or any heavy items.

Safety Instructions

- Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
- Do not operate this product if it has a damaged cord or plug.
- Do not use this product if the power socket is loose.
- Do not carry this product by its supply cord or use the cord as a handle.
- When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.

Attention


- 
- 
- Do not fall asleep while massager is on.
 - Do not use this product while you are drunk or physically uncomfortable.
 - Don't stand on the product or put a heavy weight on it.
 - Do not use this product for more than 30 minutes.
 - Do not use this product if the cloth cover is worn out.
 - Do not use this product within an hour before or after a meal.
 - When you are using the rollover function, when the flip angle has reached its maximum, please do not continue to flip hard, or it may damage the product.
 - When you are using the product at a low room temperature, do not suddenly use the heat function. You can adjust the room temperature to normal room temperature, before starting the heat function.
 - When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
 - Do not use this product continuously if you find any components damaged or the product in an abnormal condition.
 - When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
 - Stop using the product if there is thunder and lightning.
 - Do not allow children to play with the appliance.
 - Cleaning or maintenance shall not be made by children without supervision.
 - The appliance have heated surface. Person insensitive to heat must careful when using the appliance.

Warranty Claims

1. Warranty is not transferable.
 2. This warranty is effective only if the product is purchased from OGAWA.
 3. Proof of purchase (original receipt) is required for all warranty repairs.
 4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
 5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
 6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
 7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
 8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
 9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.
- 
- 

Product Parts

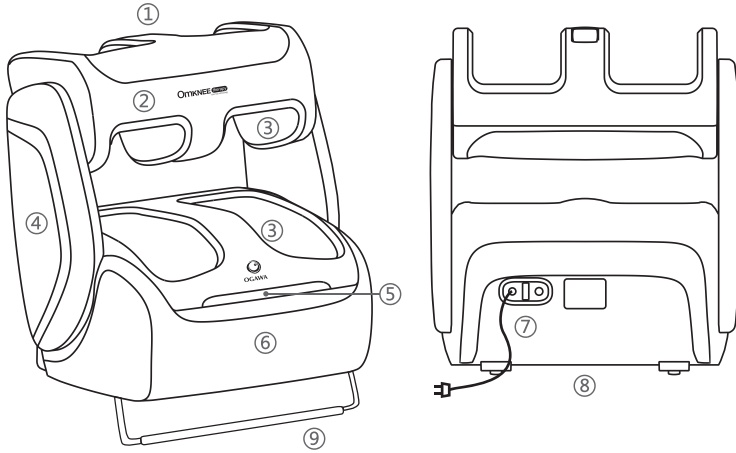
- 1 Touch screen panel
- 2 Rollover-style leg massage area
- 3 Detachable fabric cover
- 4 Side cover
- 5 Decorative plate
- 6 Feet massage area
- 7 Power switch
- 8 Two sets of transportation wheels on the ground
- 9 Stand

Product Features And Operation Illustrations


- Specially designed for knee care; U-style wrapped air compressing massage with a comfortable warm heating function, to help relieve knee pains and discomfort.
- Rotate it as you wish, – this is the first developed innovative rotation structure with 225° rotation space/11 adjustment levels, enabling you to enjoy it freely.
- New programmes: Tone/Revive/OmKNEE/Reflexology programmes are designed for different groups, to protect your legs, relieve your pressure, and ease your body and mind.
- Personalized lifting adjustment design, and 46mm automatic lifting function to satisfy different-height groups.

Product operation illustrations

Leg part

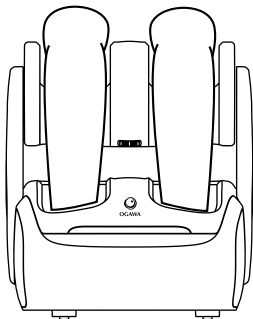Knee part

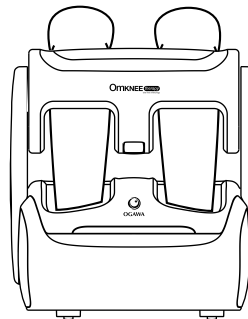
Remarks: If you would like to adjust a different angle in the process of massage, please make sure you have first turned off all the functions, before moving out your legs and then doing the necessary adjustment.


Product Operation Guide

Illustration for touch screen control panel

Leg part

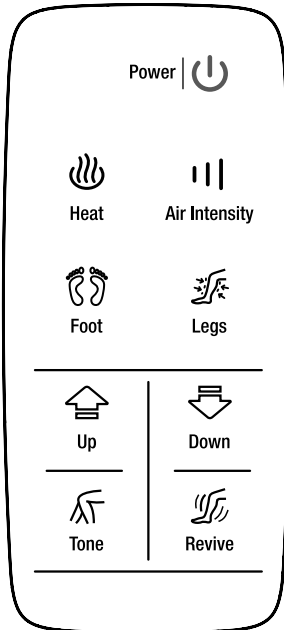Knee part

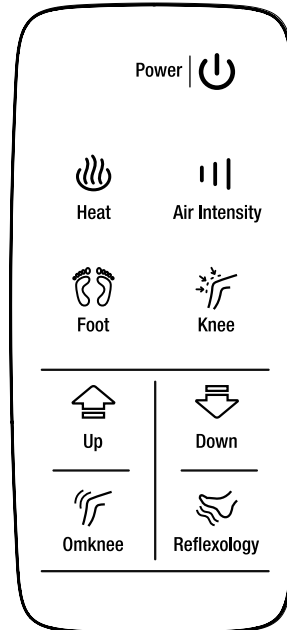① Power switch

Turn the power supply on or off, the yellow light should start to blink once the power is switched on; press the power key once again. When the blue light is on, all the other key functions will start to work.

② Intensity key Air Intensity

Control the intensity of the massage. When the product starts working, we can adjust the intensity accordingly (the default level is medium)

Intensity illustration: green light means weak / blue light means medium / red light means strong.

Product Operation Guide

③ Keys for Parts Option

- “Feet”  Air pressure and sole rolling massage for feet; enable the foot massage with the blue indicator light on;
Foot
- “Legs”  Air pressure massage for legs; enable the leg massage with the blue indicator light on;
Legs
- “Knees”  Air pressure massage for knees; enable the leg massage with the blue indicator light on;
Knee

④ Up/Down Keys

- “Up”  Turn up the leg massage region; a short buzz for the highest point;
Up
- “Down”  Turn down the leg massage region; a short buzz for the lowest point;
Down

⑤ Mode Keys

- “Tone”  Air pressure massage for legs; when enabled, the blue indicator light is on;
Tone
- “Revive”  Multi-site massage; when enabled, the blue indicator light is on;
Revive
- “OmKNEE”  Heating air pressure massage for knees; when enabled, the blue indicator light is on;
OmKnee
- “Reflexology”  Kneading massage for soles; when enabled, the blue indicator light is on;
Reflexology

⑥ Heating Key

Control the heating function. The heating effect is controlled separately; when enabled, the red or blue indicator light is on;

Product Operation Guide

a) Instructions before use


Please store this product in a dry place, do not expose the product to direct sunlight. Put the carpet under the product to minimize the working noise.

b) How to use it

This product uses only the AC 220V-240V, 50/60Hz power supply.

Make sure your indoor voltage fits this specification parameter.



- Connect the AC power cord to the socket.
- Turn on the power switch at the back of product to "I", when the red light is on, the yellow light on control panel blinks, the product is in standby mode.

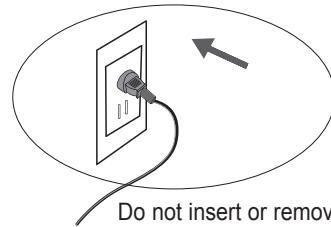Do not insert or remove the plug with a wet hand

Before plugging in the socket, make sure turn the main power switch to "0" 

- Choose your desired massage area to begin the massage.

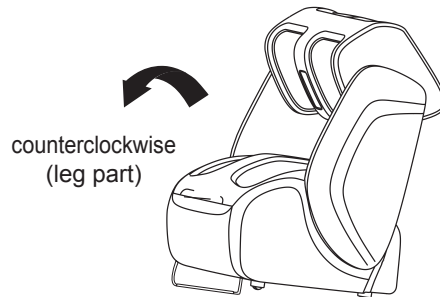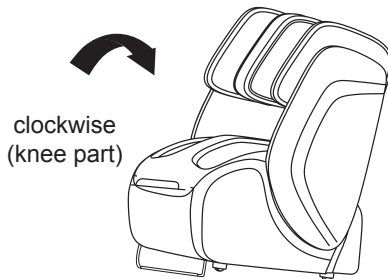(Remarks: Product rollover angle is limited to 11 levels of adjustment, and a 225 degree rotation. Once it reaches the peak point, please turn it gently and continuously, to avoid product damage.)

Product Operation Guide

- Place your legs in the product, before pressing the switch key on screen. When the blue light is on, the other functions will follow. According to the indicators on the control panel, you can choose different massage modes.

Area selection button : Choose the corresponding massage area key. When the blue light is on, three options can be selected; Press once again, to turn off the function.

Auto Program buttons : Choose the corresponding mode key. When the blue light is on, the default heating function starts to work. When the blue light is on, it is at medium level. The heat function can also be controlled. Three options can be selected.

Up/Down button : Pressing Up or Down key to adjust the angle you need. When the angle reaches its peak or lowest point, the buzzer will sound to warn you.

- When you start the massage, the product will be countdown in 15minutes. After 15 minutes, the massage function is turned off, and it will recover to standby mode.
- When the product is on standby, press the main switch “0” to cut off the power.
- Transportation wheels

There are transportation wheels installed at the base of the product. This is for moving and storage. Please hold the two sides of the product tightly, press down to lift up front portion before moving the product around.

- This appliance is design to use in sitting condition only.


Product Maintenance

a) Storage

- Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- Do not twist the power plug.
- Do not hang the power wire and pull the product.
- Avoid sharp items from destroying the surface of the product.
- Please don't store the machine under direct sunlight or in places of high temperature.

b) Clean

- Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- Use a damp cloth to clean the product. Ensure the product is completely dry before reuse.
- The detachable fabric cloth can be washed in a normal washing way.

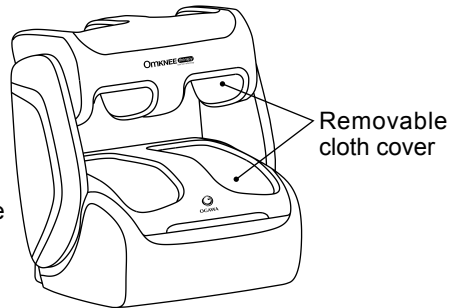*Attention

- Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

c) Maintenance

- Keep the product clean and store in a dry environment.
- When the product is not in use for a long time, please put a dustproof cover on the product.

Troubleshooting

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Conditions	Possible reason	Solutions
It doesn't work	<ol style="list-style-type: none"> 1. Plug isn't connected correctly to the socket. 2. The power switch at the back is off. 3. Switch key on control panel isn't on. 4. The fuse blows out. 	<ol style="list-style-type: none"> 1. Push the plug firmly into the socket. 2. Switch on the power switch at the back of the product. 3. Switch on the on/off key on the panel. 4. Insert a tooth pick in the small hole of "fuse" under the switch power, before taking out the fuse.
Stop working suddenly	<ol style="list-style-type: none"> 1. The power switch is off. 2. Massage time is over 15 minutes. 3. If it is working for a long time, the motor will be protected against overheating. 4. If you put heavy pressure on the rollers and the rollers do not work well. 	<ol style="list-style-type: none"> 1. Check that the power plug is connected well. 2. Restart the power supply switch. 3. After the product is cooled down, please restart the ON/OFF power. 4. Switch off the power supply, and turn it on again.

Product Specifications

Product name	OmKNEE Therapy, Total Foot Reflexology
Model number of the product	OF 2002A
Voltage	AC 220V-240V,5A , 50/60Hz
Power consumption	115W
Timer	15 minutes

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F, YHC Tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel: +85 2 3465 6262
Fax: +85 2 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM
OGAWA Vietnam
171 Nguyen Thi Thap Street
Tan Phu Ward
Dist 7
Ho Chi Minh City, Vietnam
Tel : +84 8 3775 3093 / 3775 3094
Fax: +84 8 3775 3095

INDONESIA
PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / + 632 844 6851
Fax: +632 812 0421

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