User Manual OF2005 使用手册

All information is correct at the time of print

Safety Instructions	03
Warranty Claims	06
Product Parts	07
Product Features And Operation Illustrations	80
Product Operation Guide	09
Product Maintenance	13
Troubleshooting	14
Product Specifications	15
Contact Us	16

Thank you for purchasing OGAWA OmKnee 3.0 leg massager. Please read the Safety Precautions carefully before using the leg massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the degree of danger. This book should be kept available for future reference after reading. In case of user manual get lost, user manual can also be downloaded from our website.

Safety Instructions

a) Important safety instructions

- To ensure a safe and proper operation, please read the instruction manual carefully before using this product.
- Use this product only as described in this manual.
- Always unplug this product from the electrical outlet immediately before cleaning.
- Do not use this product if you have sensitive skin, neuropathy
 Don't let children or people with mental disabilities use the heating function without proper supervision. Improper use may be dangerous.
- Do not operate this product under a blanket or when covered. Over-heating can occur and cause fire, electric shock or injury.
- Do not cover this product in order to avoid electrical faults or damages due to excessively high motor temperature.
- Stop the massage immediately if you feel any discomfort.
- Never operate this product if it has been dropped, has a damaged cord or plug;
 experiences a malfunction, or if it has been dropped in water. Please contact an authorized service centre to do the necessary check and repair.
- Do not use the components that are not recommended by OGAWA.
- Do not insert any objects into any openings in the product.
- If you want to cut off the power supply, please make sure you have turned off all the
 operating controls, as well as remove the power plug from the socket.
- Please do not use this product as a medical device. If you need a professional diagnosis or treatment, please consult a professional medical person.

b) Operating Environment Precautions

- Do not use this massager in humid or extremely dusty environments to avoid malfunction.
- Do not use this massager at a temperature above 40 degree centigrade.
- Do not use this massager outdoors.
- Avoid placing or operating this massager in direct sunlight or near a heating appliance, stove.

c) User precautions

- Make sure that this product is far away from any hot surfaces.
- Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, spinal scoliosis, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc.
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.

Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

d) Precautions for power supply

- Do not turn on the switch when your hands are wet.
- Do not immerse in or spray water or liquid on the control power switch.

e) Cord/Plug Precautions

- Do not connect or remove the plug with wet hands. Otherwise you will suffer an
 electric shock or a product malfunction.
- Do not pull the wires when you pull out the power plug. Please hold the plug itself.
- As the default operating time is 15 minutes, all the functions will stop working after 15 minutes. Turn off the main ON/OFF switch when not in use
- Keep the cord away from hot surfaces.
- Do not place the power plug under the product or any heavy items.

Safety Instructions

- Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
- Do not operate this product if it has a damaged cord or plug.
- Do not use this product if the power socket is loose.
- Do not carry this product by its supply cord or use the cord as a handle.
- When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.

Attention

- Do not fall asleep during massage.
- Do not use this product while you are drunk or physically uncomfortable.
- Don't stand on the product or put a heavy weight on it.
- Do not use this product for more than 30 minutes.
- Do not use this product if the cloth cover is worn out.
- Do not use this product within an hour before or after a meal.
- When you are using the rollover function, when the flip angle has reached its maximum, please do not continue to flip hard, or it may damage the product.
- When you are using the product at a low room temperature, do not suddenly use
 the heat function. You can adjust the room temperature to normal room temperature,
 before starting the heat function.
- When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
- Do not use this product if you find any components damaged or the product in an abnormal condition.
- When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
- Stop using the product if there is thunder and lightning.

Warranty Claims

- 1. Warranty does not cover accessories and upholstery.
- 2. Warranty is not transferable.
- 3. This warranty is effective only if the product is purchased from OGAWA.
- 4. Proof of purchase (original receipt) is required for all warranty repairs.
- 5. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
- 6. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
- Warranty does not cover accessories and attachments which do not belong to the massage unit.
- 8. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
- 9. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
- If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.


The instructions state that

- the appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction
- 2. children being supervised not to play with the appliance
- For a part of class III construction supplied from a detachable power supply unit, the instructions state that the appliance is only to be used with the unit provided

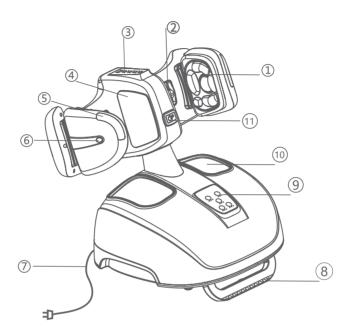- ① Massage swing arm
- ② Leg massage head
- 3 Legs operation button
- 4 Leg air bag
- ⑤ Ambiance light
- 6 Battery active key

- 7 Plug
- (8) Stand
- 9 Foot operation button
- 10 Foot fabric cover
- 11 Detach key

Detach key


Hold this button inward with one hand, and lift the upper section up with the other hand at the same time, so that the upper and lower sections can be separated.

Battery activate key

A. When using the leg section alone, you need to press this key to activate the battery, before press the power button to start the massage function.

B, This key is invalid when the leg massager is connected to foot massager.

Product operation illustrations

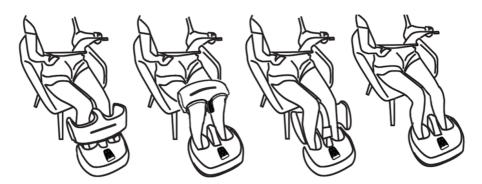Remarks: If you would like to adjust a different angle in the process of massage, please make sure you have first turned off all the functions, before moving out your legs and then doing the necessary adjustment.

Product show are for illustration purpose only

Product Operation Guide

- O Power / Mode
- In combination successful mode After the power is turned on, the "Power" light on the lower section will flashes and then turn off, and all the buttons on the lower section are invalid
- Or when the upper battery is used (separated), press this button to turn on the machine, the 15 minute countdown will begin

P1 (defaullt) blue LED
P2 green LED
P3 red LED

press again to turn off the machine in the cycle. Note: In any operating mode, long press this button to turn off all functions.

In the startup state, press this key to adjust the temperature: Low > Medium > High; Press again to turn off heat function, the indicator is off.

In the startup state, press this key to adjust the pressure intensity: Low > Medium > High; The indicator is off.

િએ Knead

Under the startup state, press this key to adjust the kneading intensity: Low > Medium > High; Press again to turn off air function, the indicator is off.

When upper and lower part is connected, only upper button is valid. When both part is separated, both button can be use

Display window

(battery level, mode program, countdown)

A. In standby mode, the battery level is displayed. When upper part power is low, it only can be charge when the combination is sucessful.

- B. In activating, when switching mode, the mode program P1, P2, P3 is displayed
- C. During massage in progress, it display 15 minute countdown

Separate condition, the indicator light flashes after the power is turned on, press this button to turn on the machine and start the M1 mode program (default); Press it again to turn off the machine. Note: In any operating mode, press this button to turn off all functions.

When the machine is on, press this key to adjust the temperature: Low > Medium > High; press again to turn off.

When the machine is on, press this key to adjust the pressure intensity: Low > Medium > High; Press again to turn off.

🚱 Knead

When the machine is on, press this key to adjust the speed: Low > Medium > High; The indicator is off.

(B) Mode

When Product is activated, press this button to change massage mode

M1 (default) blue LED
M2 green LED
M3 red LED

press again to turn off the machine in the cycle

Product Operation Guide

a) Fasten instructions:

When both part conjointed, please do the follow steps to check if they are well connected:

1:when connected, the upper sense light turns to green, and you will hear a prompt tone.

2:When the upper and lower buckle position is fastened, you will hear a click tone, indicating that the upper part has cooperated with the buckle position.

3:Lift the upper section gently to check whether it is loose from the buckle. If it is loose, please operate again.

b) Instructions before use


Please store this product in a dry place, do not expose the product to direct sunlight. Put the carpet under the product to minimize the working noise.

c) How to use it

This product uses only the 220-240V $^{\sim}$, 50-60Hz power supply.

Make sure your indoor voltage fits this specification parameter.

• Connect the AC power cord to the socket.

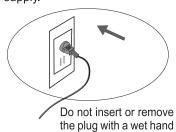- Place your legs in the product, before pressing the switch key on screen. When
 the blue light is on, the other functions will follow. According to the indicators on the
 control panel, you can choose different massage modes.
- When you start the massage, the product will be countdown in 15minutes. After
 15 minutes, the massage function is turned off, and it will recover to standby mode.

Product Maintenance

a) Storage

- Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- Do not twist the power plug.
- Do not hang the power wire and pull the product.
- Avoid sharp items from destroying the surface of the product.
- Please don't store the machine under direct sunlight or in places of high temperature.

b) Clean

- Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- Use a damp cloth to clean the product.
 Ensure the product is completely dry before reuse.
- The detachable fabric cloth can be washed in a normal washing way.

*Attention

 Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

c) Maintenance

- Keep the product clean and store in a dry environment.
- When the product is not in use for a long time, please put a dustproof cover on the product.

Troubleshooting

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Conditions	Possible reason	Solutions
It doesn't work	 Plug isn't connected correctly to the socket. There was no power in the upper section The panel key is not open 	 Push the plug firmly into the socket. The upper part can only be used after charging Open the panel
Stop working suddenly	 The power switch is off. Massage time is over 15 minutes. If it is working for a	 Check that the power plug is connected well. Restart the power supply switch. After the product is cooled down, please restart the ON/ OFF power. Switch off the power supply, and turn it on again.

Product Specifications

Product name	Omknee3.0
Model number of the product	OF2005
Voltage	220-240V~,50-60Hz
Power consumption	60W
Timer	15 minutes

•

Contact Us

MALAYSIA

Healthy World Lifestyle Sdn Bhd No. 22 Jalan Anggerik Mokara 31/47 Kota Kemuning, 40460 Shah Alam Selangor Darul Ehsan, Malaysia

Tel: +1800-88-4688 Fax: +603-5121 4386

SINGAPORE

OGAWA Health-Care Pte Ltd 163, Kallang Way #07-11 SINGAPORE 349256

Tel: +065-6846 8588 Fax: +065-6846 7677

INDONESIA PT Healthy World Lodan Center Block O No 20 Jalan Lodan Raya No 2 Jakarta Utara 14430

Tel: +62 21-6983 1123 Fax: +62 21-6983 1124

HONG KONG, CHINA

OGAWA Health Care International(HK) Ltd 31/F,YHC tower,1Sheung Yuet Road

Kowloon Bay Kowloon

Hong Kong

Tel:+852 3465 6262 Fax:+852 2410 0119

VIETNAM

OGAWA Vietnam Co Ltd 2nd floor,Petroland Tower, 12 Tan Trao, Tan Phu ward District 7

Ho Chi Minh

Tel: +848 54133 222/54133 223

Fax: +848 5416 1736

PHILIPPINES

O-Healthcare Solution Phil. Inc. 2303, 4th Floor, G&A Building, Chino Roces Avenue Extension, Makati City, Philippines

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OGAWA Vietnam Co Ltd 2nd floor,Petroland Tower, 12 Tan Trao, Tan Phu ward District 7 Ho Chi Minh