# *mobile*seat **iBACK 2**

SMART MASSAGE CUSHION

Ogawa is grateful  
to have received  
these awards:

PUTRA  
BRAND AWARDS 2010  
SILVER  
In association with  
MALAYSIA'S MOST VALUABLE BRANDS

Ogawa Wellness Hubs: ·Malaysia·Singapore·Vietnam·Philippines·Indonesia  
·Myanmar·India·China·Hong Kong·Taiwan·Turkey·Pakistan  
·Saudi Arabia·Mauritius·Russia·France·Canada

User Manual 使用手册  
OG 1001

# CONTENTS

---

Safety Instructions	2
Warranty Claims	4
Main Unit & Accessories	5
Controller Diagram	6
Product Operation Guide	7
Product Features	9
Product Maintenance	10
Troubleshooting	11
Product Specifications	12
Contact Us	13

# Safety Instructions

---

## i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Place the appliance in a safe location.
- Do not use or sit on the appliance when driving the car.
- Do not expose the appliance to a heater or direct sunlight.

## ii. Unsuitable user

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease
- Pregnancy
- Any health concerns

## iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Do not use this product before sleep.
- If the fabric cover is torn.
- If liquid is spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

## iv. Caution and warning

- Unplug the appliance from the outlet when not in use.

# Safety Instructions

---

- Do not use any attachments that are not recommended by OGAWA.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.
- Do not insert your hands or any object into the massager's vents and gaps.
- Use the heated surfaces carefully. It should not be used by children or incapacitated persons. May be dangerous.
- Do not lift or carry this appliance by the power cord or any cord.
- This product is only intended for home use only. Commercial use is prohibited.
- Do not cover the appliance during operation.
- Do not crush appliance and avoid it from sharp folds.

## v. Note

- Please read these safety notes thoroughly beforehand to ensure the safety and correct use of the appliance.
- The appliance is not intended for use by persons suffering from physical, sensory or mental disabilities, or persons who lack experience and knowledge about the product. They need to be supervised and instructed by a person responsible for their safety.
- Extensive usage could lead to excessive heating and shorter product lifespan. Should this occur, discontinue usage and allow the unit to cool before operating.
- Some levels of noise during usage are normal due to vibration.
- Consult the doctor before use if the person using it is under medication or with a medical condition. Please do not use the product if you are not feeling well.
- Do not increase the room temperature abruptly if the appliance is used in the cold room. It is advisable that the temperature should increase gradually to a normal level.
- If the appliance is stored in a cold place or in a warm environment, please wait for one hour before use to reach the optimum massaging level.

## Safety Instructions

---

- Do not remove the cover of the appliance. Touching the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power cord from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.

## Warranty Claims


---

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other nonresidential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.


## Main Unit & Accessories

---

- 1 Back cushion
- 2 Controller
- 3 Adaptor
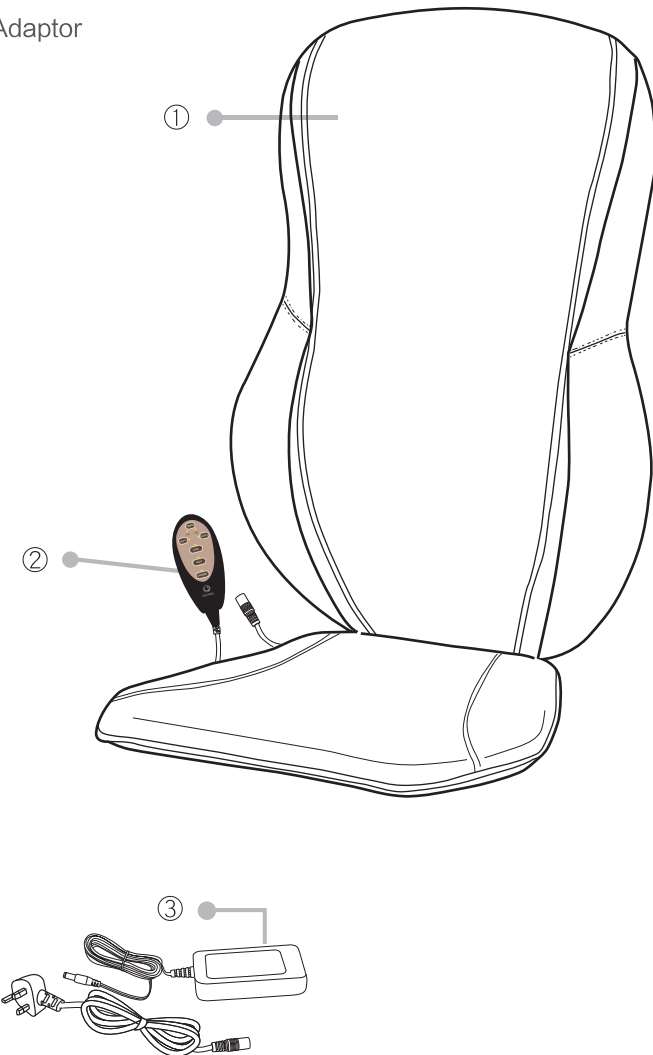
## Controller Diagram

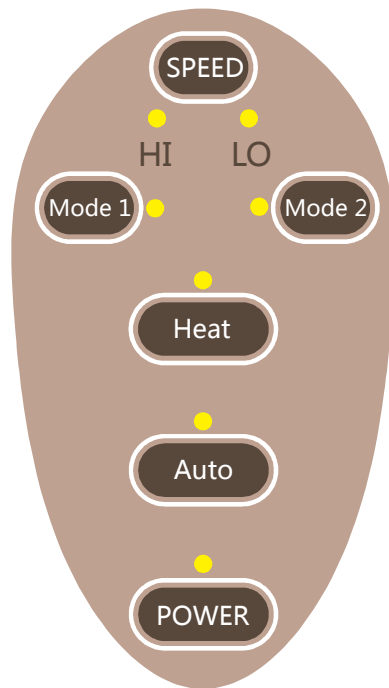**NOTE: There is a 15-minute auto shut off on the unit for your safety.**

## Product Operation Guide

### i. Preparation before use

- A screw is applied on the back of the product to protect the massager apparatus during shipment. Remove the screw before use.
- Use the equipped wrench to remove the screw and insert the power cord provided into the opening.

### ii. Placing the massager

- The massager is equipped with back strap that allow user to faster onto most chairs.
- Place the massager on a chair, wrap the two ends of the strap round the back of the chair and secure with Velcro.

## Operation instruction

### i. [Power]

- Insert the power plug into the electric outlet. Massage roller will go back to standby position.
- Press POWER to start massage in low speed, timer is preset at 15 minutes.
- Press POWER to stop.

# Product Operation Guide

---

## ii. [AUTO]

Press AUTO to begin auto program.

## iii. [Mode1]

Press Mode 1, the rollers will move in same rotation.

## iv. [Mode2]

Press Mode 2, the rollers will move in different rotation.

## v. [Speed]

Press SPEED to adjust the intensity of massage.

## vi.[Heat]

Press HEAT to start the heat function.

# Product Features

---

1. 3 massage zones for your selection. Concentrate on your upper back, lower back or your entire back.
2. Heat function to increase the massage effects and improve blood circulation  
Light and portable design.
3. Easy to use remote controller.
4. Auto shut-off after 15 minutes.
5. LED light indicators.

## Effects by using the appliance

1. Relieves the muscle ache and pain.
2. Relieves the tension and reduces the stress.
3. Improves blood circulation.

When not in use, press the [Power] button, and please don't cut off the power immediately until the rollers have resumed to the bottom and the indicator light is off.

## After use

Press the 'Power' button after use and the massager will stop operating. Massage rollers will resume to their original position. Unplug the appliance from the power outlet.

## Product Maintenance

---

### To store

- Keep the massager in the box or in a safe, dry and cool place.
- Keep the massager away from sharp edges or pointy objects that might cut or puncture the product's fabric surface.
- To avoid damage, do not wrap the power cord around the unit.
- Do not hang the unit with the controller cord.

### To clean

- Avoid immersing unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Keep the appliance away from water or any liquid.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish or paint thinner to clean the product.

## CAUTION

- It is recommended for every massage not more than 15 minutes.
- Turn off the power switch immediately if user feels any discomfort.
- Unplug the appliance from the power outlet when not in use.
- If out of power, please turn off the power switch to avoid the massager from continue operating once the electricity is restored.
- This product is auto shut-off after 15 minutes for safety purpose.

## Troubleshooting

---

his massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

Conditions	Solutions
Not operating.	<ul style="list-style-type: none"><li>• Use the product when electricity is restored.</li><li>• Reconnect the power code.</li><li>• Turn on the power switch.</li></ul>
Sudden stop.	<ul style="list-style-type: none"><li>• Use the product when electricity is restored.</li><li>• Reconnect the power cord.</li></ul>

## Product Specifications

---

Product Name	iBACK2
Model number of the Product	OG1001
Power source	220-240V~
Frequency	50-60Hz
Voltage	24V 2.5A
Power consumption	60W
Dimensions (Length/width/height)	L608 X W170 X H434 (mm)
Weight	4.8KG

## Contact Us

---

**MALAYSIA**  
 Healthy World Lifestyle Sdn Bhd  
 No.22 Jalan Anggerik Mokara 31/47  
 Kota Kemuning, 40460 Shah Alam  
 Selangor Darul Ehsan, Malaysia  
 Tel: +603-5121 4286  
 Fax: +603-5121 4386

**HONGKONG**  
 OGAWA Health Care  
 International (HK) Ltd  
 31/F, YHC Tower,  
 1 Sheung Yuet Road, Kowloon Bay,  
 KLN, Hong Kong S.A.R.  
 Tel : +85 2 2428 8823  
 Fax: +85 2 2410 0119

**SINGAPORE**  
 OGAWA Healthy-Care Pte Ltd  
 27 Tampines Industrial, Avenue 5  
 T5 @ Tampines  
 Singapore  
 Tel : +065-6846 8588  
 Fax: +065-6846 7677

**VIETNAM**  
 OGAWA Vietnam  
 2nd Floor , Petroland Tower ,  
 12 Tan Trao , Tan Phu Ward ,  
 District 7 , Ho Chi Minh  
 Tel : +84 8 3775 3093 / 3775 3094  
 Fax: +84 8 3775 3095

**INDONESIA**  
 PT Healthy World  
 Lodan Center Block O No 20  
 Jalan Lodan Raya No 2  
 Jakarta Utara 14430  
 Tel : +62 21-6983 1123  
 Fax: +62 21-6983 1124

**PHILIPPINES**  
 O-Healthcare Solution Phil. Inc.  
 2303, 4th Floor, G&A Building,  
 Chino Roces Avenue Extension,  
 Makati City, Philippines  
 Tel: +632 889 0676 / +632 844 6851  
 Fax: +632 812 0421

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