
# LENI TOUCH

Handheld Vibration Massager

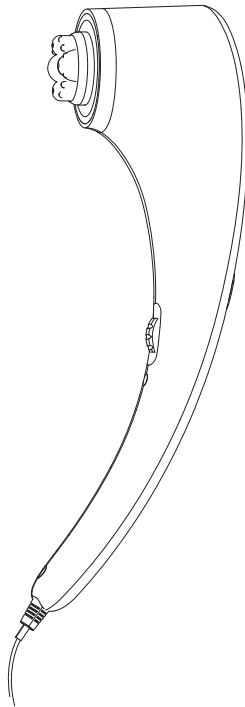User Manual  
OL 0361 使用手册

# LENI TOUCH

Handheld Vibration Massager

User Manual

All information is correct at the time of print

# CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Controller Diagram	06
Product Operation Guide	07
Product Maintenance	08
Troubleshooting	09
Product Specifications	09
Contact Us	10

# Safety Instructions

## i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not use the appliance on the bed.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink. Place the appliance in a safe location.
- Do not sit on the appliance.
- Do not expose the appliance to a heater or direct sunlight.

## ii. Health precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease
- Pregnancy
- Any health concerns

## iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Do not use this product before sleep.
- If the cover is break.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

## iv. Caution and warning

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by OGAWA.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.

# Safety Instructions

- Do not drop or insert any object into any opening.
- Do not lift or carry this appliance by the power cord or any cord.
- This appliance is designed only for personal use, and not for professional purposes.
- Do not cover the appliance during operation.
- Do not crush the appliance and avoid any sharp folds.

## v. Note

- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- The appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. They need to be supervised and instructed about its usage by a person responsible for their safety.
- Extensive usage could lead to excessive heating and a shorter product lifespan. Should this occur, discontinue the use and allow the unit to cool before operating again.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.


# Safety Instructions

- The appliance is only to be used with power supply provided.
- The appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.

## Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

# Main Unit & Accessories

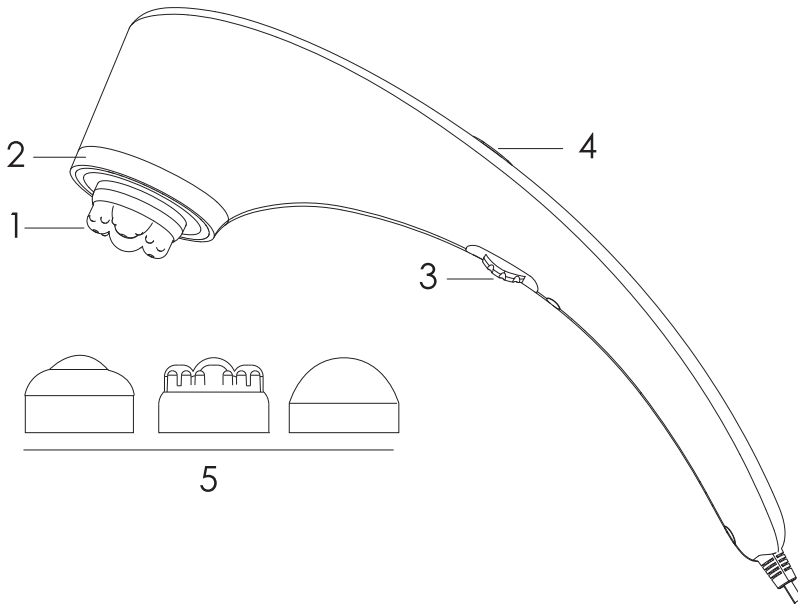1. Massage Nodes
2. Heating Areas
3. On-Off Switch/Speed Control  
(Infinitely Adjustable Massage Intensity)
4. Heating Function Switch
5. Interchangeable Massage Attachments.

# Product Operation Guide

## **i. Preparation before use**

1. Before Using The Appliance
  - Select the proper massage head.
  - Plug the appliance to the power outlet.

## **ii. Operation guide**

1. To turn on, twist up the speed control until a "click" sound is heard.
2. After the machine been activated. The tapping speed can be change by twist up the speed control.
3. To turn off, twist down the speed control until a "click" sound is heard. The machine will off immediately.
4. Press the heating button to turn on the heat, the heating area will light up showing the heat is activate. Press again, to turn off

## **iii. Benefits from using the appliance**

1. Relieves muscle aches and pain.
2. Relieves tension and reduce stress.
3. Improves blood circulation.

## **iv. After Use**

Unplug the appliance from the power outlet.

# Product Maintenance

## TO STORE

- Place the massager in this box or in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the surface.
- Please wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

## TO CLEAN

- Avoid immersing the unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to go inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish or paint thinner to clean the product.
- Do not attempt to repair the product. For servicing, please contact OGAWA.

## CAUTION

- Avoid tearing on the appliance surface.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.

# Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

Problem	Cause	Remedy
Massage head does not move.	Appliance not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Device does not move.	Switch on the appliance.
	Overheating protection tripped.	Disconnect device from mains power. wait at least 15 minutes. plug in again and switch on.

# Product Specifications

Product Name	Leni Touch
Model number of the product	OL0361
Adaptor voltage	220-240V~
Power frequency	50-60Hz
Power consumption	20W
Dimensions (Length/width/height )	L425mm*W93mm*H140mm
Weight	1 kg±0.1 kg

# Contact Us

**MALAYSIA**  
Healthy World Lifestyle Sdn Bhd  
No. 22 Jalan Anggerik Mokara 31/47  
Kota Kemuning, 40460 Shah Alam  
Selangor Darul Ehsan, Malaysia  
Tel : +603-5121 4286  
Fax: +603-5121 4386

**HONG KONG, CHINA**  
OGAWA Health Care International (HK) Ltd  
31/F, YHC Tower, 1 Sheung Yuet Road  
Kowloon Bay  
Kowloon  
Hong Kong  
Tel : +852 3465 6262  
Fax: +852 2410 0119

**SINGAPORE**  
OGAWA Health-Care Pte Ltd  
27 Tampines Industrial Avenues S,  
T5 @Tampines  
Singapore 528623  
Tel : +065-6846 8588  
Fax: +065-6846 7677

**VIETNAM**  
OGAWA Vietnam Co Ltd  
2nd floor, Petroland Tower,  
12 Tan Trao, Tan Phu ward  
District 7  
Ho Chi Minh  
Tel: +84 8 3775 3093/3775 3094  
Fax: +84 8 3775 3095

**PHILIPPINES**  
O-Healthcare Solution Phil. Inc.  
2303, 4th Floor, G&A Building,  
Chino Roces Avenue Extension,  
Makati City, Philippines  
Tel : +632 889 0676 / +632 844 6851  
Fax: +632 812 0421

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OGAWA Health Care International (HK) Ltd  
31/F,YHC Tower,1 Sheung Yuet Road  
Kowloon Bay  
Kowloon  
Hong Kong  
Tel: +852 3465 6262  
Fax: +852 2410 0

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OGAWA Health-Care Pte Ltd  
27 Tampines Industrial Avenue 5,  
T5 @ Tampines  
Singapore 528623  
Tel: +065-6846 8588  
Fax: +065-6846 7677

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OGAWA Vietnam Co Ltd  
2nd floor,Petroland Tower,  
12 Tan Trao,Tan Phu ward  
District 7  
Ho Chi Minh  
Tel: +84 8 3775 3093/3775 3094  
Fax: +84 8 3775 3095

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O-Healthcare Solution Phil. Inc.  
2303, 4th Floor, G&A Building,  
Chino Roces Avenue Extension,  
Makati City, Philippines  
Tel : +632 889 0676 / +632 844 6851  
Fax: +632 812 0421

Ogawa is grateful  
to have received  
these awards:

Ogawa Wellness Hubs : Malaysia . Singapore . Vietnam . Philippines . Indonesia . Myanmar . India . China  
Hong Kong . Taiwan . Turkey . Pakistan . Saudi Arabia . Mauritius . Russia . France  
Canada