
WonderEase+

Portable Neck and Back Massager

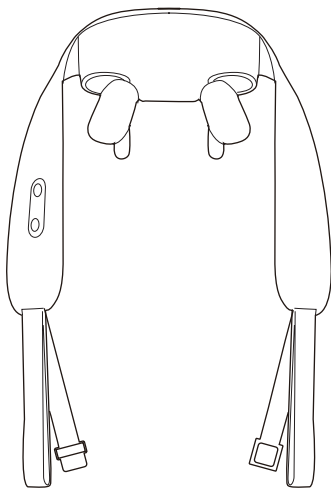WonderEase+

Portable Neck and Back Massager

User Manual

All information is correct at the time of printing

Contents

Important Safeguards	2
Instruction for Use	5
Warranty Claim	7
Product Specification	8
Operating Instruction	10
Contact Us	12

Important Safeguards

IN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

DANGER- TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug the appliance from the electrical outlet immediately after complete charged.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.

WARNING - TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE OR INJURY TO PERSONS:

- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by OGAWA; specifically any attachments not provided with the unit.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to OGAWA Service Center for examination and repair.
- Keep cord away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.

Important Safeguards

- Do not use the product when charging
- To disconnect, turn all controls to the “off” position, then remove plug from outlet.
- DO NOT use outdoors.
- NEVER operate the appliance with the air openings blocked. Keep the air openings free of lint, and hair.
- NEVER operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

SAVE THESE INSTRUCTIONS

Caution—Please read all instructions carefully before operating.

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if
 - You have a pacemaker
 - You are pregnant
 - You have any concerns regarding your health
- Not recommended for use by Diabetics.
- Never leave the appliance unattended, especially if children are present.
- Never cover the appliance when it is in operation.
- Do not use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.

Important Safeguards

- Never use this product directly on swollen or inflamed areas or skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- The product is equipped with overheating protection circuit, if the power is cut-off, please wait until the product is cooled down.
- DO NOT use this product before bed. The massage has a stimulating effect and can delay sleep.
- NEVER use this product while in bed
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the owner half of their body.
- This unit should not be used by children or invalids without adult supervision.


Instruction for Use

Charging and energy reminder function:

1. When charging with an external switch power supply, the product cannot be turned on and the charging power supply must be unplugged to start the massage.
2. During the charging process, the charging indicator light breathes red.
3. When fully charged, the charging indicator light stays green.
4. When the battery energy is $<10\%$ or the NTC fails or the motor is blocked and shut down, the red indicator light flashes, and after 5 times, the product shuts down.

*Stabbing shutdown: If more than two stalls are detected within 5 seconds, the product shuts down.

*Note: If the massager is continuously charged for more than 6 hours and the red light still breathes (in charging state), please unplug the power adapter to avoid danger.

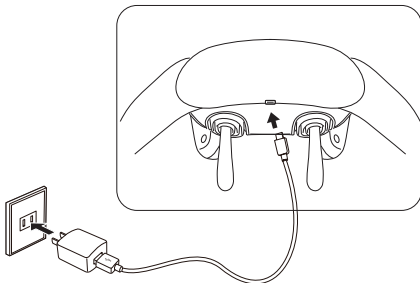Instruction for Use

Light status	function	short press
Orange light always on	kneading + heat (L)	<p>*Mode button: Every time you switch mode, the indicator light flashes once, and then the corresponding color indicator light stays on.</p> <p>*Heating button: Every time you switch mode, the indicator light flashes twice, and then the corresponding color indicator light stays on.</p>
Red light always on	kneading +heat (H)	
Green light always on	kneading	
Red light breathing	charging	
Green light always on	full charged	
Red light flash	low battery, NTC error	

Warranty Claim

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA or any authorized distributors.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repairs are not covered by this warranty.
6. Warranty does not cover accessories and add-ons, which do not belong to the product.
7. Warranty is not effective to rental, business, commercial, institutional, or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitutions in lieu of repair or replacement.
10. This appliance is not intended for use by persons(including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.Children should be supervised to ensure that they do not play with the appliance.
11. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Product Specification

Product name	WonderEase+
Model no.	OL0843
Charging voltage	5V
Charging current	2A
Power rating	20W
Operating voltage	7.4V
Battery capacity	2500mAh
Charging time	est 3.5hrs

This product is for home use only

*When using this product, you must use a switch power cord certified by our company.

*This product is only equipped with a Type-C charging cable for charging the product. You need to purchase a charging head by yourself. The charging head is recommended to use 5V 2A specifications.

*If the product is not used for a long time, please charge the product within half a year to activate the battery.

Correct usage method Preparation before use

* Please check whether the product and accessories are complete and whether there is any damage.


* Please massage gently when starting to use.

* If there is any abnormality in the body during use, please stop using it.


* The use time of this product is set to 10 minutes at a time, and it will stop automatically after the time is up. It is recommended not to exceed 30 molecules per use.

Product mix


Product Body

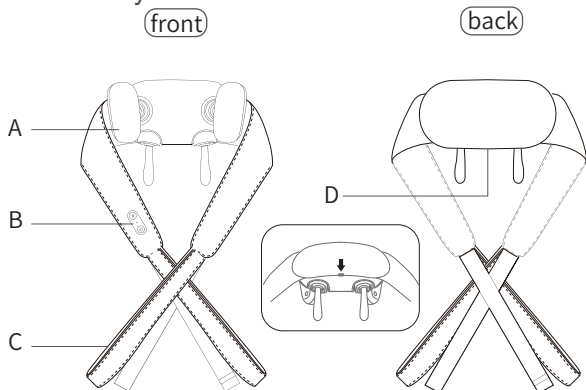- A. Stainless steel heater
- B. Controller
- C. Buckle
- D. Type-C connector (bottom)

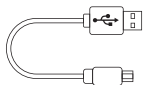Type-C cable

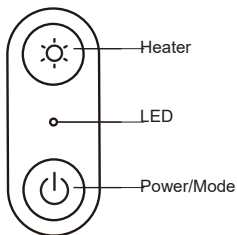Operating Instruction

Power/Mode Button:

- Power on: Press and hold this button for 2 seconds to power on, default kneading mode 1, heating at low temperature 40 degrees, orange light always on, timer for 10 minutes.
- Power off: Press and hold this button for 2 seconds to power off, or timer for 10 minutes to automatically power off, massage head reset to open state.
- Mode selection: Power on default mode 1, short press this button to switch modes; the switching mode sequence is: mode 2 - mode 1 - mode 2, and so on. If only the mode is turned on without heating, the green light is always on.

Mode 1: kneading.

Mode 2: vibration + kneading.

Heating button: (Low temperature 40℃, High temperature 45℃)

- In the power-on state, default low temperature, the indicator light is orange; press this button to enter high-grade heating, the indicator light is red, and then press it to turn off the heating.

Switching order: low heating-high heating-turn off heating, and so on.

Low heating: the orange indicator light is always on.

High heating: the red indicator light is always on.

Heating off: the green indicator light is always on.

Operating Instruction

To Store

- Place massager in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, do not wrap the power cord around the unit. Do not hang the unit by the hand control cord.

To Clean

Unplug the unit and allow it to cool before cleaning. Spot clean only with a soft, slightly damp sponge. Never allow water or any other liquids to come into contact with the unit.

- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner to clean.
- Do not attempt to repair the Neck Massager.

Contact Us

MALAYSIA

Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +1800-88-4688
Fax: +603-5121 4386

HONG KONG, CHINA

OGAWA Health Care International (HK) Ltd
31/F, YHCTower, 1 SheungYuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel :+852 3465 6262
Fax:+852 2410 0119

SINGAPORE

OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM

OGAWA Vietnam
Co Ltd PETROLAND TOWER
2nd floor, 12 Tan Trao Street,
Tan Phu Ward, Dist 07,
HCMC, Vietnam
Tel:+84 28 54133 222/223

PHILIPPINES

O-Healthcare Solution Phil.Inc
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676/ +632 844 6851
Fax: +632 812 0421

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