HANDSON

NECK MASSAGER

User Manual

All information is correct at the time of printing

Contents

Important Safeguards	2
Instruction for Use	5
Warranty Claim	6
Product Specification	7
Functional Instruction	8
Maintenance	9
Contact Us	10

Important Safeguards

IN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

DANGER- TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug the appliance from the electrical outlet immediately after complete charged.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.

WARNING - TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE OR INJURY TO PERSONS:

- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by OGAWA; specifically any attachments not provided with the unit.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to OGAWA Service Center for examination and repair.
- Keep cord away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.

Important Safeguards

- Do not use the product when charging
- To disconnect, turn all controls to the “off” position, then remove plug from outlet.
- DO NOT use outdoors.
- NEVER operate the appliance with the air openings blocked. Keep the air openings free of lint, and hair.
- NEVER operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

SAVE THESE INSTRUCTIONS

Caution—Please read all instructions carefully before operating.

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if
 - You have a pacemaker
 - You are pregnant
 - You have any concerns regarding your health
- Not recommended for use by Diabetics.
- Never leave the appliance unattended, especially if children are present.
- Never cover the appliance when it is in operation.
- Do not use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.

Important Safeguards

- Never use this product directly on swollen or inflamed areas or skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- The product is equipped with overheating protection circuit, if the power is cut-off, please wait until the product is cooled down.
- DO NOT use this product before bed. The massage has a stimulating effect and can delay sleep.
- NEVER use this product while in bed
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the owner half of their body.
- This unit should not be used by children or invalids without adult supervision.

Instruction for Use

Charging Instruction

LED color	LED condition	Status
White	Always on	battery capacity >25%
Orange	Flash	Low, please charge
Orange	Breathe	charging
Green	Always on	Full charge

* This product cannot be used while charging. Please unplug the power adapter after charging is completed.

* The product uses fast charging technology. It is normal for the host and power adapter to heat up slightly during charging.

* Please use the charging cable provided by our company to charge the product and pay attention to electricity safety.

* Check the power supply. If the power supply voltage is out of the specified range, this product cannot be charged.

Warranty Claim



1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA or any authorized distributors.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repairs are not covered by this warranty.
6. Warranty does not cover accessories and add-ons, which do not belong to the product.
7. Warranty is not effective to rental, business, commercial, institutional, or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitutions in lieu of repair or replacement.
10. This appliance is not intended for use by persons(including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.Children should be supervised to ensure that they do not play with the appliance.
11. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Product Specification


Product Name	Hands-on neck massager
Model	OL0844
Working Voltage	7.4V ===
Power Consumption	15W
Battery Capacity	2000mAh
Weight	0.9kg
Size	205x213x93mm

Remark : This is rechargeable product, build with type C connection

Functional Instruction

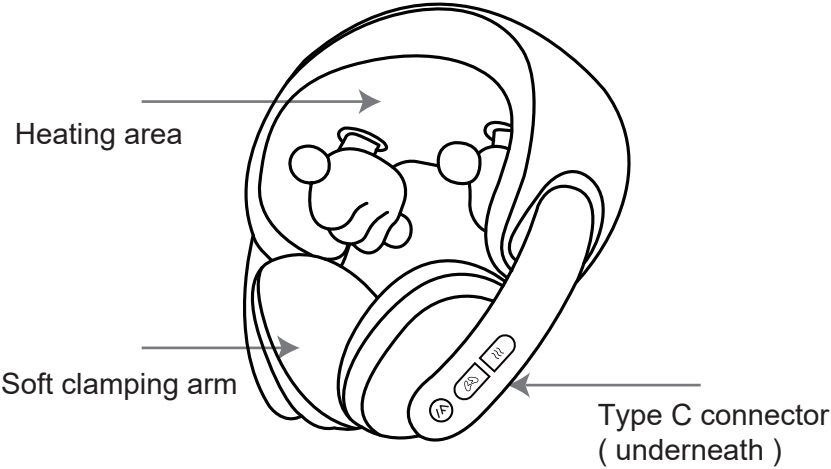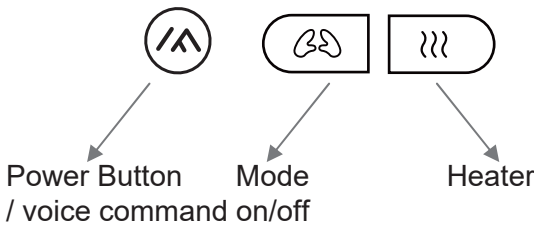USB cable

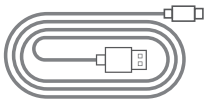Maintenance

Operating Instruction

【 Power 】: Press and hold 1.5 sec to turn on, press and hold 1.5 sec to turn off

【 Voice prompt 】: Double-click to open. Double-click again to close.

【 Mode 】: 5 massage mode available, short press to change
energetic > relaxing > neck care
> wake up > kneading

【 Heating 】: 2 Mode available, short press to change
Low > High > Turn Off

To Store

- Place massager in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, do not wrap the power cord around the unit. Do not hang the unit by the hand control cord.

To Clean

Unplug the unit and allow it to cool before cleaning. Spot clean only with a soft, slightly damp sponge. Never allow water or any other liquids to come into contact with the unit.

- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner to clean.
- Do not attempt to repair the Neck Massager.

Contact Us

MALAYSIA

Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : 1800 88 4688

HONG KONG, CHINA

OGAWA Health Care International (HK) Ltd
31/F, YHCTower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel : +852 3465 6262
Fax: +852 2410 0119

SINGAPORE

OGAWA Health-Care Pte Ltd
163, Kallang Way #07-11
SINGAPORE 349256
Tel : +065-6846 8588

VIETNAM

OGAWA Vietnam
171, Nguyen Thi Thap Street
Tan Phu Ward
Dist 7
Ho Chi Minh City, Vietnam
Tel : +84 8 3775 3093 / 3094
Fax : +84 8 3775 3095

PHILIPPINES

O-Healthcare Solution Phil.Inc
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676/ +632 844 6851
Fax: +632 812 0421

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