

# Comfi Back

Lower Back Air Traction massager




## User Manual OM 1003

### Safety Cautions


- To ensure the safety and proper use of the product, please read the instruction carefully before use. Carefully examine the covering before each use.
- Use this massager only for its intended use as described in this manual.
- For adult use only.
- Consult your physician before using this massager.
- Always unplug this massager from the electrical outlet immediately after use and before cleaning.
- Do not operate under blanket or pillow. Excessive heating can occur and cause electric shock or injury to persons.
- To avoid the motor temperature is too high to cause trouble or damage, please do not cover the product.
- Stop operation immediately if you experience any discomfort.
- Never operate this unit if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. If the unit's power supply cord or plug is damaged, do not attempt to fix it by yourself.
- Do not use attachments not recommended by manufacturer.
- Never drop or insert an object into any opening.
- The product is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- Do not use products in damp or dusty environment so as to avoid product malfunction.
- Do not place the product outdoors.
- Following person, please consult with the doctor before use.
  - Who uses a medical electrical devices, which is embedded in the body, such as pacemakers.
  - Who suffers from malignant tumors
  - Who suffers from heart disease.
  - Who is pregnant, after delivery or in a period monthly illness
  - Who has a high degree of peripheral obstruction caused by diabetes and intuitive damage
  - Who is suffering from skin trauma
  - Whose body temperature over 38°C (fever) such as: acute inflammation symptoms (burnout, aversion to cold, blood pressure changes)
  - Who is in serious or weak period
  - Who has osteoporosis, vertebral fractures, sprains, muscle strain, acute disease
  - Whose spine with abnormal or incliner.
  - Who are treated in a medical institution.
- Do not use wet hands to connect or disconnect the power plug from the electrical outlet.
- Do not place, drop or submerge in water or other liquid.
- Do not operate for more than 30 minutes continuously.
- Never leave the unit unattended, especially if children are present.
- Do not fall asleep during product operation.

### Product Description

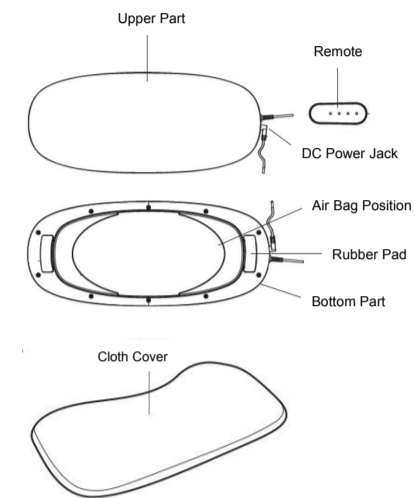

### How to Operate

Connect the adapter and DC power Jack of the unit, then connect the power supplier.

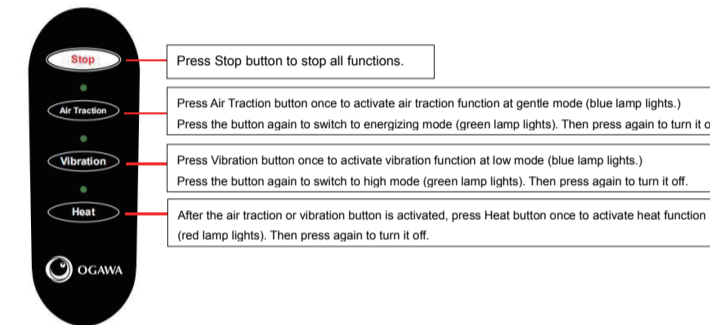

Note: The product will automatically power off in 10 minutes.

### Maintenance

#### Cleaning

- Before cleaning, turn the massager off and unplug the plug from the power outlet.
- Wipe dust and dirt off the product with a dry, soft cloth after use.
- Do not wash or soak product.

#### Storage

- Unplug the plug from the power outlet.
- Place the massager in a safe, cool and dry place when not in use.
- Do not put a heavy load on the unit.
- Store out of the reach of children.

### Troubleshooting

This page summarizes the most common problems you could encounter with the product. If you are unable to solve the problem with the information below, please contact us immediately. Please do not disassemble or repair the product yourself.

Problem	Possible Cause	Action
The product does not work.	The product is not connected with the adapter properly.	Check the connection between the adapter cord plug and the unit jack.
	The operation button is off.	Turn the operation button on.
	The product has just finished its 10 minutes auto timer session.	If you want another session, please press the operation button to turn the product on.
	The product has been used continuously for a long time, then the overheat protection has been activated.	Unplug the product and Let it cool down before using it again.
The product is abnormally warm.	An excessive load has been applied on the product.	Reduce the excessive load and start again.
	The product is likely to be overused.	Unplug the product and Let it cool down before using it again.
Noise can be heard from the product.	Noise is made by the mechanical parts of the product.	This is normal and does not mean that the product is faulty.

### Specification

Name of Product:	Comfi Back
Model No.:	OM 1003
Rated Voltage:	AC 100-240V / 50/60Hz
Power Output:	DC 12V / 1.2A
Power Consumption:	14.4W
Auto Timer:	10 minutes
Product Dimension:	Approx. 400*172*57mm
Product Weight:	Approx. 1.06kg (adapter excluded)

### Contact us

MALAYSIA  
Healthy World Lifestyle Sdn Bnd  
No. 22 Jalan Anggerik Mokara 31/47  
Kota Kemuning, 40460 Shah Alam Selagor Darul Ehsan, Malaysia  
Tel: +603-5121 4286  
Fax: +603-5121 4386