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

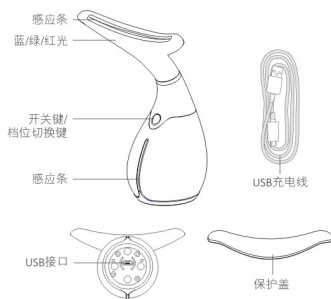
HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F, YHC tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel: +852 3465 6262
Fax: +852 2410 0119

VIETNAM
OGAWA Vietnam Co. Ltd
2nd floor, Petroland Tower,
12 Tan Trao, Tan Phu ward
District 7
Ho Chi Minh
Tel: +84 8 3 775 3093/3775 3094
Fax: +84 8 3775 3095

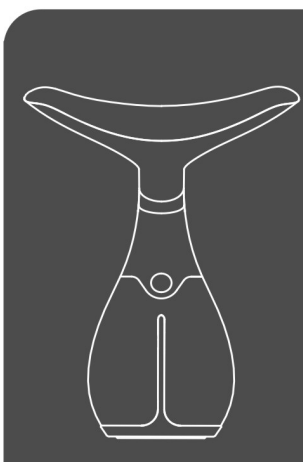OY 0160
INSTRUCTION MANUAL


STRUCTURE:

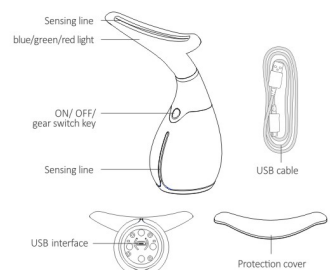FEATURES:

- Newly adopted the bionics design concepts, shaped in dolphin tail, looks novelty and classic.
- The massage head is specially designed for shoulder, neck and other parts of body, streamline outlook fit for ergonomic theory, can help transporting the energy to every part of skin.
- The massage heads heating in 45℃ to activate the skin, matching with three IPL acoustic mode in low, medium and high speeds, and blue, green, red 3 color light massage modes, can promote the blood circulation, stimulate the neck acupuncture points, clear the meridians, improve circulation and promote metabolism, promote the discharge of long-accumulated toxic substances & waste, slow skin cell aging, to dilute the neck profile, less swelling, relax neck, make acupuncture treatments, rapidly eliminate fatigue, restore the role of physical fitness.
- Three sonic modes as low speed (blue), medium speed (green), high speed (red), which are all made by intelligent technology, after selecting mode, when massage head contacts the skin, the vibration starts.
- With mode memory function and 60 seconds intimate reminder function, it can easily adjust the massage time.
- USB cable charging, safety and environmental protection.
- Mini & portable, can be used anywhere any time.

INSTRUCTION OF THE 3MODES

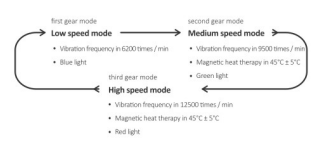Blue light is the coldest light, which can promote the synthesis of protein and collagen, can activate the skin, tighten the loose skin and shrink the pores, is good for oily skin and sensitive skin.

Green light can improve the oxygen alternate use function of the cell, promote the skin surface microcirculation, and dredge the lymphatic & edema, meanwhile it is with the function of calming skin, balancing skin and keeping stability of the skin.

Red light can increase the cell vitality, accelerate the blood circulation, and promote the growth of fibroblasts and collagen, has significant effect in smoothing fine lines, wrinkles, firming skin, color dodging, and fading the scar.

APPLICABLE USER

This Massager is suitable for all persons.

HOW TO USE:

1. Apply the skin care product to the neck.
2. Select the mode according to personal need, hold the massage head lift neck from bottom.
3. Back neck also need massage, use device massage on accupunctual points by ways of pressing, squeezing, stroking, pushing etc.

HOW TO CLEAN

1. Long press ON/OFF for 1 seconds to turn off the product.
2. Wipe the product with cleansing cloth, do not use the detergents.
3. Close the protection cover, and put the product on a dry and ventilated place.

CAUTIONS

1. Please read the product manual before using, and keep the manual for future use.
2. USB cable charging, charging from low battery station to saturation needs 2 hours, for first gear and second gear, fully charged allows 3 hours of continuous use.
3. When battery is in shortage, flash in red light, when it is under battery charging, flash in blue light, when blue light is up on indicates charge in saturation.
4. The product has mode memory function, when switch it on, it will enter the mode you chose last time automatically.
5. The product will be turned off automatically if no use for more than 2 minutes. (If there is any operation within 2 minutes, it will remain standby status).
6. It is not allowed to use this product when it is under charging.
7. The product has the function of heating, approximate temperature is in 45℃. Different people are in some difference when sensing the temperature, some are likely to feel hot. After you apply the essence oil, massage cream on skin, the product will automatically adjust the temperature to most comfortable.
8. Please stop to use the product if you feel uncomfortable or skin complaint and consult a doctor.
9. The product is not waterproof; please do not place it in water.
10. Put the product on a dry and ventilated place, to avoid direct sunlight, and make sure no corrosive chemicals around.
11. The cleansing cloth using high-quality fleece material, can be washed repeatedly.
12. Keep out of reach of children, to avoid any danger.

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn. Bhd
No. 22 Jalan Anggerik Mokara 31/47,
Kota Kemuning, 40460 Shah Alam,
Selangor Darul Ehsan, Malaysia
Tel: +603-5121 4286
Fax: +603-5121 4386

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5@Tampines
Singapore 528623
Tel: +065-6846 8588
Fax: +065-6846 7677

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F, YHC tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel: +852 3465 6262
Fax: +852 2410 0119

VIETNAM
OGAWA Vietnam Co. Ltd
2nd floor, Petroland Tower,
12 Tan Trao, Tan Phu ward
District 7
Ho Chi Minh
Tel: +84 8 3 775 3093/3775 3094
Fax: +84 8 3775 3095