

封面底面


封面材质:120g双胶纸、内页材质:80g双胶纸 封面颜色:金色PANTOEN 871C 内页颜色:金色PANTOEN 871C+PANTONE COOL GRAY 8C 装订方式:平装胶装





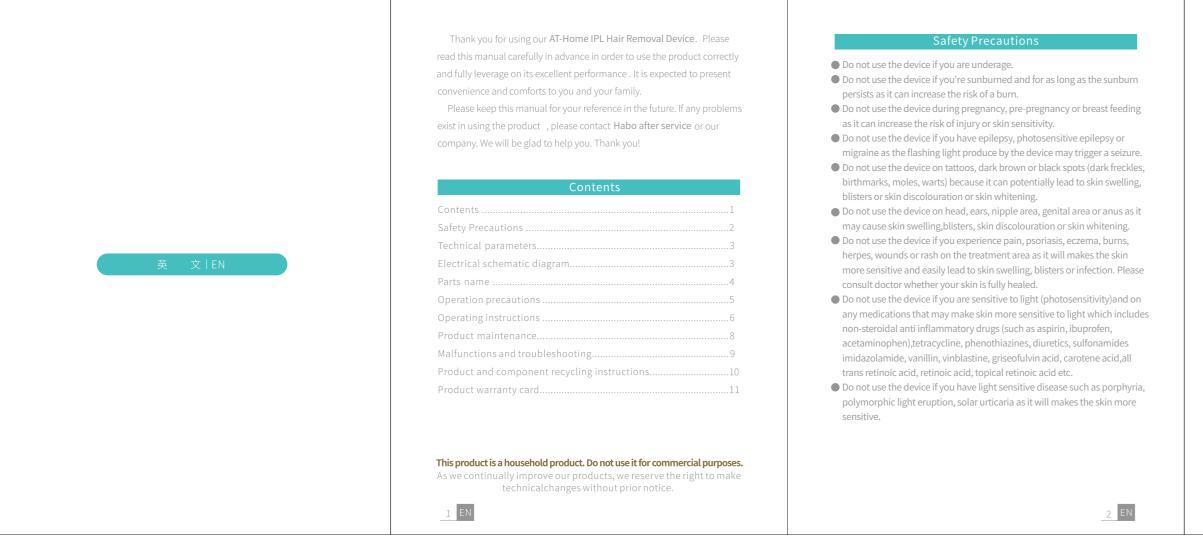

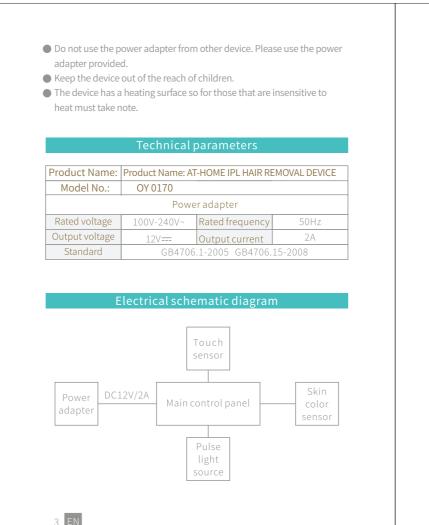

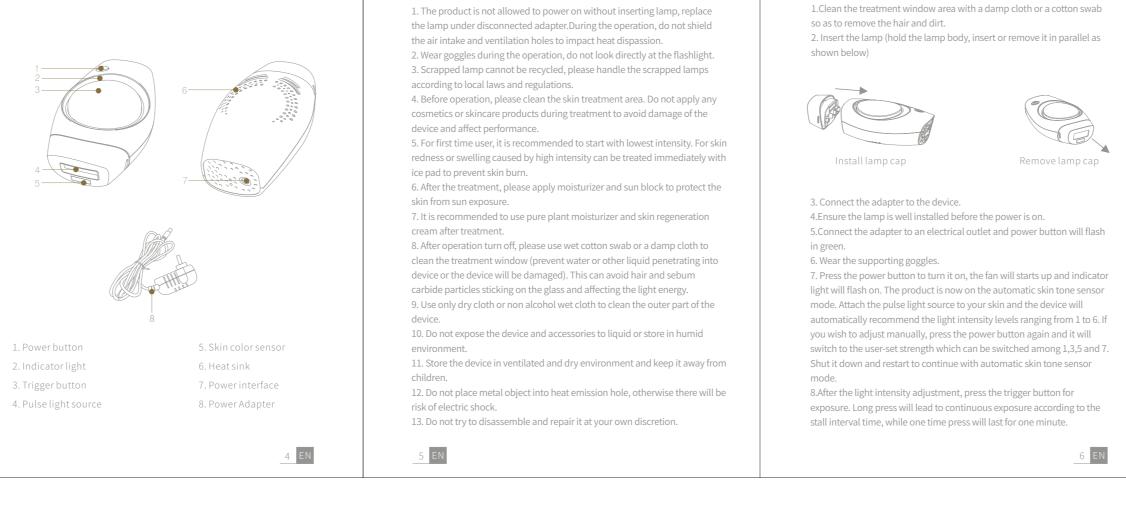

the air intake and ventilation holes to impact heat dispassion. 2. Insert the lamp (hold the lamp body, insert or remove it in parallel as 2. Wear goggles during the operation, do not look directly at the flashlight. 3. Scrapped lamp cannot be recycled, please handle the scrapped lamps cosmetics or skincare products during treatment to avoid damage of the device and affect performance.

* first time user, it is recommended to start with lowest intensity. For skin redness or swelling caused by high intensity can be treated immediately with Install lamp cap ice pad to prevent skin burn. 6. After the treatment, please apply moisturizer and sun block to protect the skin from sun exposure. 3. Connect the adapter to the device. 4.Ensure the lamp is well installed before the power is on. 7. It is recommended to use pure plant moisturizer and skin regeneration cream after treatment. 5. Connect the adapter to an electrical outlet and power button will flash 8. After operation turn off, please use wet cotton swab or a damp cloth to clean the treatment window (prevent water or other liquid penetrating into 6. Wear the supporting goggles. device or the device will be damaged). This can avoid hair and sebum 7. Press the power button to turn it on, the fan will starts up and indicator carbide particles sticking on the glass and affecting the light energy. light will flash on. The product is now on the automatic skin tone sensor 9. Use only dry cloth or non alcohol wet cloth to clean the outer part of the mode. Attach the pulse light source to your skin and the device will automatically recommend the light intensity levels ranging from 1 to 6. If 10. Do not expose the device and accessories to liquid or store in humid you wish to adjust manually, press the power button again and it will environment. switch to the user-set strength which can be switched among 1,3,5 and 7. 11. Store the device in ventilated and dry Shut it down and restart to continue with automatic skin tone sensor 8. After the light intensity adjustment, press the trigger button for 12. Do not place metal object into heat emission hole, otherwise there will be

exposure. Long press will lead to continuous exposure according to the

stall interval time, while one time press will last for one minute.

9.The green light on power button flashes if it fails to detect the skin tone so make sure the lamp is attached firmly on the skin. During operation, if the green light flash rapidly, it is a warning sign. IPL isn't suitable for the darkest skin tones. 10.The device comes with 7 intensity which includes LED 1, LED 2, LED 3, LED 4, LED 5, LED 6, LED 7. 11. Long press of power button or if the device is in idle status for more than 3 minutes, it will shuts down, and the indicator light and fan will turn off as well and the power button will be flashing slowly.

RECOMMENDED IPL TREATMENT						
1st-2nd week	3rd-4th week	4th week & above				
3 times a week	2 times a week	1 time a week				
lair begins to fall	Hair grow is much	Any regrowth can				

by 75% noticeable touch-up

* Results may vary from individual to individual

device or accessories. Storage 1.Unplug the device from the power outlet and clean the device. 2.Keep the device in a cool, dry place away from sunlight. 3. Do not put heavy objects on the device. out and is reduced finer and less be managed with

1.Do not clean the device with water. Instead, use a dry cloth or soft cloth

with little soap to wipe. Any contact with water can damage its machinery.

2. Do not use alcohol, banana oil, gasoline, nail polish remover and other

corrosive substance to wipe the device as it may cause damage to the

 1. Light outlet failing to be attached to skin or dark

1. Attach the light outlet to the skin or change skin ears spet. skin color leadsto failure in exposure work with power indicator light flashing(green light in Failure in exposure in latter)

work

2.More than 1-minute one -time use against skin 3. Move away from skin while pressing the flash button

4.Broken modulator tube leads to malfunction

away to rest your skin and then attach the skin again to start using it and so on.

3. Loosen the button after re-attaching the skin and press the flash button before using it again.

4. Shut down and stop using

If by following the instructions mentioned above can not solve the faults,

properly please stop using the product at once.

the outer case was damaged or breakdown, the product still can not work

Failure phenomenon	Causes	Troubleshooting		Recycle indicativ	e instru	ctions		
ail to use he product	Lamp cap and skin color Re-operate according are not calibrated and access to main unit. the power is off from the dome				ct the environment, when onsumer you have the respicic waste and send to a qua			
Failure in exposure work	1. Light outlet failing to be attached to skin or dark	to the skin or change skin care spot according to the instruction. 2. Move the product away to rest your skin and then attach the skin again to start		the product or its p with relevant nation Name and co	nal regu	ulations		
	skin color leadsto failure in exposure work with power indicator light flashing(green light in former condition, green				Hazardoı			
				Part name	Pb	Hg	Cd	
				Plastic parts (ABS, PC etc.)	0	0	0	
				Power adapter	×	0	×	
	in latter)			PCB parts*	×	0	×	
	2.More than 1-minute one -time use against skin			Metal parts (interior)	0	0	×	
	3. Move away from skin while pressing the flash			Packaging (color box,manual, nameplate, etc.)	0	0	0	
	button 4.Broken modulator tube leads to malfunction			O: Indicates the homogeneous requirement x: Indicates the	us mate t specifi	rial of the	ne comp /T 2657:	
Automatic shutdown	1.System continuous non- exposure operation for 3	Restart Shut down and stop		least one ho the limit req *: Indicates the	mogene uiremer	eous ma nt specif	terial of ied in Gl	
	minutes 2. Failure in fan	using		This product is "green life" of 10 y property designed by the compar				

Warranty provisions: 1. Consumers buy for family use and the appliances enjoy one-year machine warranty since the date of purchase in normal condition; 2.If you need free maintenance for the above-mentioned f ailuresduring the warranty period, please take invoice to the

maintenance point. 3. No free maintenance can be enjoyed in case of any of the following circumstances. However, the paidmaintenancewill be offered. a . Damage caused by improper use, maintenance and storage of the product by customers; b. No valid invoice;

c. The model does not match with the product model which requires repairing or has been altered; d. Damage caused by the force majeure; 4. In addition to the maintenance issues specified in the product, the holders can enjoy other rights within the scope of the relevant

national laws and regulations.

内页