


# 100\*70mm 上下对折, 正反两面 铜版纸128g

## User Manual

Ion Cleansing & Infusing & Cooling Device

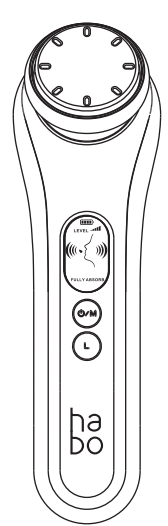

OY 0171

Dear customer, thanks for purchasing this product. Please read this user manual carefully and keep it for future reference.

**Product Introduction**

Habolon Cleansing & Infusing & Cooling Device is designed to make your daily facial clean and skin care more effective. By holding the device, it will form a current loop between itself and the human body. Features constant heating of 45°C/113 °F and high frequency vibrations, the device, based on the principle of charges repelling, opens the pores and takes out toxin/grease/dirt deep in pores. At meantime, based on the principle of opposite charges attracting, the device infuses water and nutrients into the skin, helps skin absorb them quickly, and thus improves the skin from the inside out. Besides, the device comes with a special cool head to help shrink the pores and calm the skin. Regularly using the device will bring you a clearer and more glowing face.

**Safety Warning**

For your safety, please read the user manual carefully before use and keep it for future reference. If the device is transferred to others, the user manual should be attached, too.

**Precautions**

Please observe the following precautions, otherwise it may lead to device failure or personal injury.

- Do not use the device on areas around eyes or ears and other sensitive areas.
- Do not use the device on broken skin, eczema or tumour areas.
- Do not use the device until all the metal things (such as earrings, necklace, bracelet) have been removed from the body.
- Do not use the device during charging.
- Do not keep the device in one place for more than 10 seconds.
- Do not take apart, transform or repair the device by yourself.
- Do not use facial masks other than sheet masks if you are allergic to metal.
- Keep the intensity in an acceptable range when use.
- Keep the device out of reach of children.
- Do not use the device if you are pregnant unless being directed by a physician.


**Store & Care**

- Dry the device with a clean towel after each use.
- Put the device away if not in use for a long term.
- Do not wet the device with water, wash it under faucet or immerse it in water.
- Do not use corrosive cleanser to clean the device.
- Put the device in a dry, cool shadow.
- Do not use the device if it is damaged.
- Do not use the device to strong impact.


**Package Contents**

- 1\* Ion Cleansing & Infusing & Cooling Device
- 1\* Attachment Ring
- 1\* Micro-USB Charging Cable
- 1\* User Manual

**Product Overview**

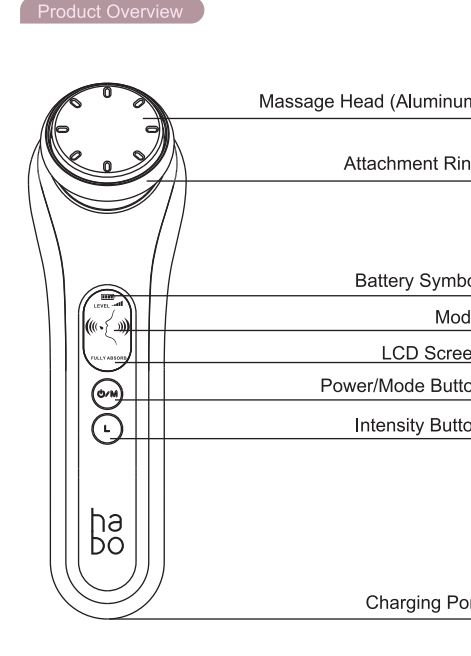

- Message Head (Aluminum)
- Attachment Ring
- Battery Symbol
- Mode
- LCD Screen
- Power/Mode Button
- Intensity Button
- Charging Port
- Cool Head (Aluminum)
- Ion Induction Area

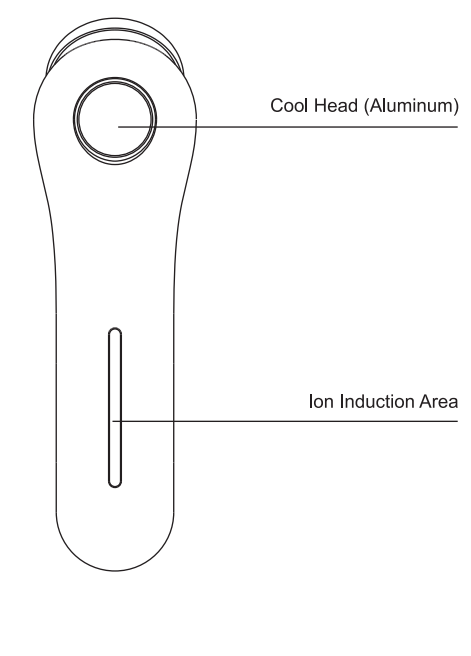

Ion Induction Area

**Product Features**

- Customized: Comes with 4 modes and 5 intensities to meet different needs
- Effective: Constant heating and high frequency vibrations for deeper clean and fuller absorption
- Special: Special cool head helps to shrink the pores, soothe and calm the skin
- Smart: Smart LCD screen shows battery life, mode and intensity clearly and directly
- Convenient: Special attachment ring makes the use of the device more convenient
- Eco-Friendly: Built-in 2000mAh Lithium battery supports long runtime, USB rechargeable
- Attentive: Low battery reminder help recharge the device in time, no interruption during use
- Ergonomic: The device body conforming to ergonomic design, feels great in hand

**Charging Methods**

- When the battery is low, the battery symbol on the LCD screen will flash, which reminds you to recharge the device in time. And the device will switch itself off in a few minutes.
- Connect the device with an adaptor (not included) via the

Micro-USB charging cable. Insert the plug into the socket. The LCD screen lights up and the battery symbol flashes, which means charging begins.

- The battery symbol will turn solid once the battery has been charged completely. Unplug the adaptor at this time.

Note -- The device doesn't work during charging.

Note -- A full charge takes about 5 hours.

Note -- To protect the battery, please recharge the device once every 3 months if not in use for a long term.

**Fit the Cotton Pad**

- Remove the attachment ring from the Massage Head.
- Put a cotton pad on the Massage Head and replace the attachment ring.

Note -- You can saturate the cotton pad with lotions before replace the attachment ring if the cotton pad is hard to be fitted.

Note -- The heat radiating from the device may be weakened if you use a thick cotton pad.

- Saturate the cotton pad with enough lotions.

Note -- The device is not waterproof. Please wipe off the remaining lotions after each use.

**Application Methods**


**Mode 1: DEEP CLEAN**

- Fit a cotton pad on the Massage Head and saturate the cotton pad with enough lotions with function of cleansing or exfoliating.
- Press and hold the Power/Mode Button for 1.5 seconds to turn the device on. The LCD screen will light up. The default mode is DEEP CLEAN mode.

The default intensity is the lowest intensity. You can press the Intensity Button to choose a certain intensity. It's recommended to start with the lowest intensity.

- Hold the device and make sure that your hand directly contacts with the Ion Induction Area.

- Gently apply the Massage Head to your face and move it slowly in a circle motion.



After use, remove the cotton pad and wipe off the remaining lotions on the Massage Head. Then wash your face with clean water.

Note -- Do not reuse a cotton pad to avoid any possible skin problems.

Note -- Do not use purified water or mineral water to replace lotions, or the facial cleansing effect will be reduced.

Note -- Ensure that your hand directly contacts with the ion induction area on the device, or the device will not work.

**Mode 2: FULLY ABSORB**

- Clean your face, and apply lotions or essence to your face.
- Press and hold the Power/Mode Button for 1.5 seconds to turn the device on. The LCD screen will light up. Press the Power/Mode Button to choose the FULLY ABSORB mode.

- The default intensity is the lowest intensity. You can press the Intensity Button to choose a certain intensity. It's recommended to start with the lowest intensity.
- Hold the device and make sure that your hand directly contacts with the Ion Induction Area.
- Gently apply the Massage Head to your face and move it slowly in a circle motion.

Note -- The mode needs no cotton pads. It's recommended to apply moisturizing skin-care products to get better effect.

Note -- When the device is hard to move, apply more lotions or essence to your face in time.

Note -- Do not keep the device in one place for more than 10 seconds. Due to the micro current, you might feel a slight tingling. It's normal and you can rest assured to use.

Note -- Ensure that your hand directly contacts with the ion induction area on the device, or the device will not work.

**Mode 3: MOISTURE**

- Apply a mask or evenly apply emulsion or cream to your face.

- Press and hold the Power/Mode Button for 1.5 seconds to turn the device on. The LCD screen will light up. Press the Power/Mode Button to choose the MOISTURE mode.
- The default intensity is the lowest intensity. You can press the Intensity Button to choose a certain intensity. It's recommended to start with the lowest intensity.
- Hold the device and make sure that your hand directly contacts with the Ion Induction Area.

- Gently apply the Massage Head to your face and move it slowly from the inside out, and bottom up.
- After use, wipe off the remaining lotions or essence on the Massage Head.

Note -- The mode needs no cotton pads or facial lotions. It's recommended to apply moisturizing skin-care products for better effect.

Note -- When the device is hard to move, apply more mask essence, emulsion or cream to your face.

Note -- Do not keep the device in one place for more than 10 seconds. Due to the micro current, you might feel a slight tingling. It's normal and you can rest assured to use.

Note -- Ensure that your hand directly contacts with the ion induction area on the device, or the device will not work.

**Mode 4: COOL**

- Press and hold the Power/Mode Button for 1.5 seconds to turn the device on. The LCD screen will light up. Press the Power/Mode Button to choose the COOL mode.
- After about 10 seconds, the temperature of the Cool Head will drop to about 10 °C/50 °F. Then apply the Cool Head to the skin and massage gently, helping shrink the pores and calm the skin.
- After use, clean the Cool Head.

Note -- This mode can work with or without skin care products. So you don't have to remove emulsion/cream/essence/mask you have applied in the previous skin-care steps.

Note -- After using the mode, the device will get hot. It's normal and you can rest assured to use.

Note -- To avoid high temperature damages the device, do not use the mode continuously. And the minimum interval between each use should be 5 minutes. If you try to reuse the mode within 5 minutes, the device will beep. It's normal.

**Specifications**

Name: Ion Cleansing & Infusing & Cooling Device  
 Material: ABS + PC + Aluminum  
 Battery Capacity: 2000mAh  
 Charging Time: 5 Hours (approx.)  
 Continuous Working Time: Cool Mode - 2 Hours (approx.)  
 Other Modes - 5 Hours (approx.)  
 Input Current: 5V/500mA  
 Output Power: MAX 11W  
 Product Size: 165\*50\*60mm  
 Product Weight: 204g  
 Package Size: 190\*87\*70mm  
 Package Weight: 415g

Contact Us:  
 MALAYSIA  
 HEALTHY WORLD LIFESTYLE SDN BHD 550407-K  
 No.22, Jalan Anggerik Makara 3/147, Kota Kemuning, 40400 Shah Alam, Selangor Darul Ehsan, Malaysia. T +603 521-4286 F +603521-4386

**Reference for Use**

	Quick Care	Simple Care
Time	4 minutes	8 minutes
Mode	MOISTURE	MOISTURE + DEEP CLEAN
Using frequency (recommend)	Once a day	Once a day
Effect	Quickly relax the skin	Deep clean the skin

	Full Care	Cool Care
Time	12 minutes	2 minutes
Mode	DEEP CLEAN + FULLY ABSORB + MOISTURE	COOL
Using frequency (recommend)	Three or four times a week	Once a day
Effect	Help the skin better absorb the nutrients for a glowing face	Help shrink the pores and calm the skin

**Troubleshooting**

Troubles	Possible Causes	Solutions
Fail to boot	Press the power button for less than 2 seconds	Press the power button for 2 seconds
	During charging	Unplug the adaptor or the device
	Low battery	Plug into the power supply to recharge
Fail to recharge (No battery symbol showed in the LCD screen)	The charging cable is not completely inserted into the adaptor	Insert the charging cable into the adaptor completely
	The charging cable is not completely inserted into the device	Insert the charging cable into the device completely
Device went hot	The charging cable is damaged	Use another charging cable
	Work continuously for too long	Turn the device off and let it rest for at least 1 hour
After using the COOL mode, the device will get hot. It's normal.	Hand has no contact with the ion induction area	Ensure hand directly contacts with the ion induction area
	Use an intensity which is too strong	Turn down to a proper intensity
Do not work	Apply not enough lotions or essence	Apply more lotions or essence
Feel a strong tingling		