EYETOUCH PLUS

WITH THERMO CARE

Ogawa is grateful
to have received
these awards:

Ogawa Wellness Hubs : Malaysia . Singapore . Vietnam . Philippines . Indonesia . Myanmar . India . China
Hong Kong . Taiwan . Turkey . Pakistan . Saudi Arabia . Mauritius . Russia . France
Canada

OY0300 User manual

EYE TOUCH **PLUS**
WITH THERMO CARE

Contact Us

Malaysia

Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel: +603-5121 4286
Fax: +603-5121 4386

Singapore

OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel: +665-6846 8588
Fax: +665-6846 7677

Indonesia

PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel: +62 21-8983 1123
Fax: +62 21-8983 1124

HongKong

OGAWA Health Care International (HK) Ltd
31/F, YHC Tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel: +852 3465 6262
Fax: +852 2410 0119

Vietnam

OGAWA Vietnam
PETROLAND TOWER
2nd Floor, 12 Tan Trao Street,
Tan Phu Ward, Dist 07,
HCMC, Vietnam
Tel: +84 8 54133 222/54133 223
Fax: +84 8 5416 1736

Philippines

O-Healthcare Solution Phil, Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel: +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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Safety instructions	3
Part name and function	5
Product structure	5
Controller	6
Product operation Guide	7
Instruction	7
Charging	8
Maintenance and Storage	9
Troubleshooting	10
Specification	11

Safety instructions

Safety 1

- ▶ Close and continuous supervision is necessary when this appliance is used by or near children, invalids or disable persons.
- ▶ If you are experiencing health complications, please consult your doctor before use.
- ▶ Do not operate this appliance if it has a damaged cord or damage charger .
- ▶ Before charging, please keep power cord, power plug and charger plug dry.
- ▶ Do not use if water is accidentally spilled onto the appliance.
- ▶ Do not use together with other medical devices.
- ▶ Do not use this product as a substitute for medical treatment.
- ▶ This product is being designed to use by one person at one time only.
- ▶ Do not use this appliance when drunk or not feeling well.
- ▶ Recommend to use this product for 15 minutes at one time.
- ▶ Do not pull the product via power cord on floor to prevent product failure.
- ▶ Those who suffer from eye disease are requested to consult a doctor prior to use.
- ▶ The following persons are requested to consult a doctor before use to reduce the risk of the user feeling sick or suffering from an injury:
 - Those who have a malignant tumor
 - Those who suffer from heart disease
 - Those who experience high blood pressure

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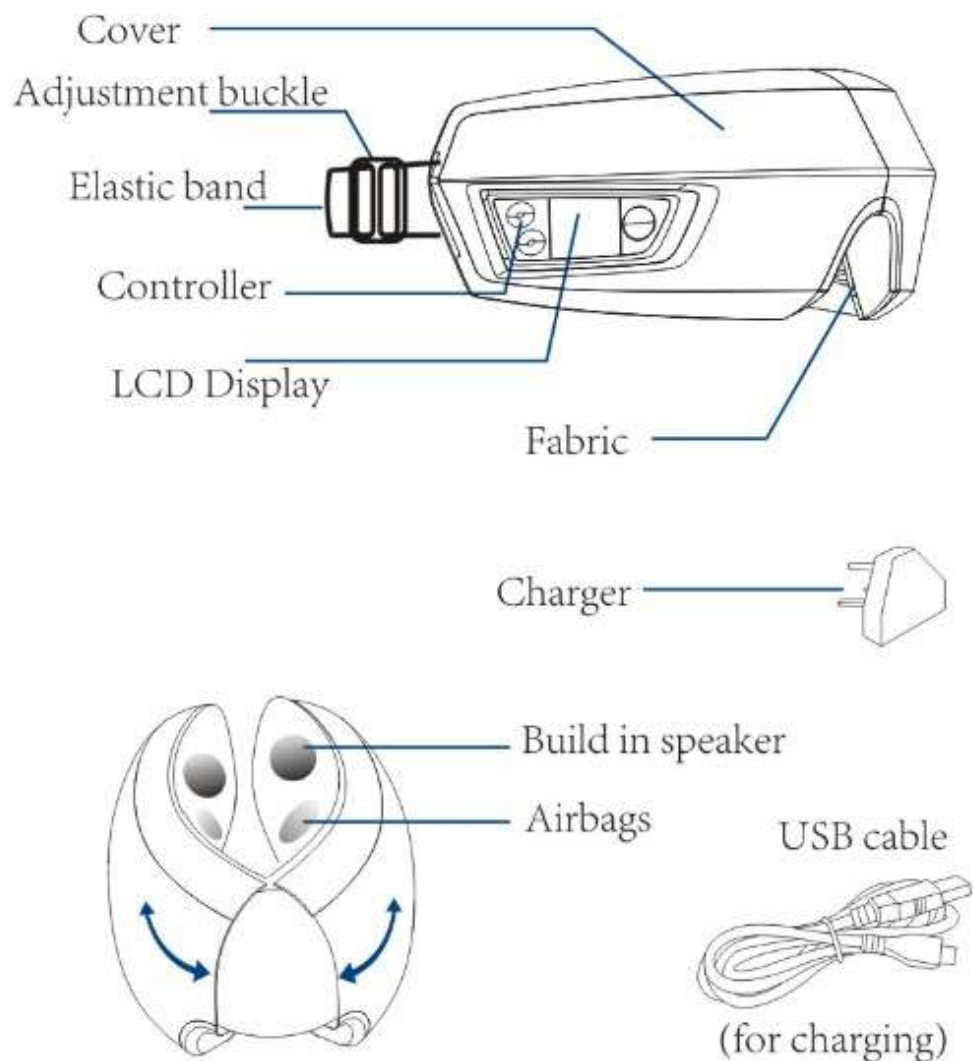
- ▶ Do not use if body temperature more than 38C (fever).
- ▶ Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction .
- ▶ All services must be done by authorized OGAWA Technicians only, unauthorized persons are not allow to dismantle or repair the appliance.
- ▶ Improper usage of electricity voltage may result in fire.
- ▶ Be sure to remove the power supply cord by pulling the power plug, not the power cord.
- ▶ Only use the charger that came together with the appliance for charging purpose.
- ▶ Take off your eye contact lenses or glasses before use this product.

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