
*mobile*seat **XE** **DUO**

DUAL-ACTION MASSAGE CUSHION

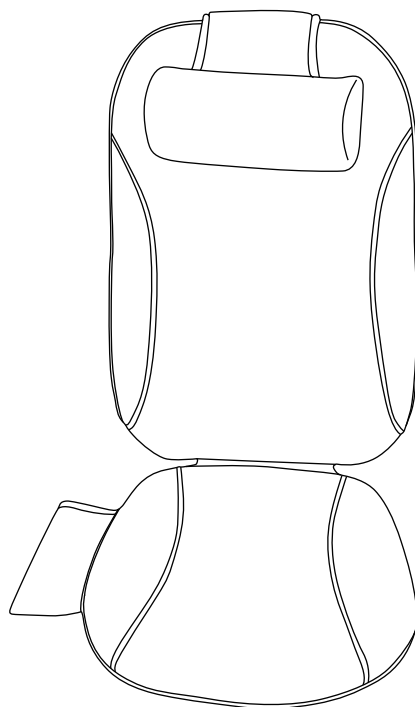User Manual
OZ 0939 使用手册

OGAWA

*mobile*seat **XE** **DUO**
DUAL-ACTION MASSAGE CUSHION

User Manual

All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	06
Main Unit & Accessories	07
Controller Diagram	08
Product Operation Guide	09
Display	12
Product Maintenance	13
Troubleshooting	14
Product Specifications	15
Contact Us	16

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Place the appliance in a safe location.
- Do not use or sit on the appliance when driving the car.
- Do not expose the appliance to a heater or direct sunlight.

ii. Unsuitable user

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease
- Pregnancy
- Any health concerns

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Do not use this product before sleep.
- If the fabric cover is torn.
- If liquid is spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

Safety Instructions

iv. Caution and warning

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by OGAWA.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.
- Do not insert your hands or any object into the massager's vents and gaps.
- Use the heated surfaces carefully. It should not be used by children or incapacitated persons. May be dangerous.
- Do not lift or carry this appliance by the power cord or any cord.
- This product is only intended for home use only.
Commercial use is prohibited.
- Do not cover the appliance during operation.
- Do not crush appliance and avoid it from sharp folds.

Safety Instructions

v. Note




- Please read these safety notes thoroughly beforehand to ensure the safety and correct use of the appliance.
- The appliance is not intended for use by persons suffering from physical, sensory or mental disabilities, or persons who lack experience and knowledge about the product. They need to be supervised and instructed by a person responsible for their safety.
- Extensive usage could lead to excessive heating and shorter product lifespan. Should this occur, discontinue usage and allow the unit to cool before operating.
- Some levels of noise during usage are normal due to vibration.
- Consult the doctor before use if the person using it is under medication or with a medical condition. Please do not use the product if you are not feeling well.
- Do not increase the room temperature abruptly if the appliance is used in the cold room. It is advisable that the temperature should increase gradually to a normal level.
- If the appliance is stored in a cold place or in a warm environment, please wait for one hour before use to reach the optimum massaging level.
- Do not remove the cover of the appliance. Touching the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power cord from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.

Warranty Claims


1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

- ① Back cushion
- ② Controller
- ③ Car adaptor
- ④ Adaptor

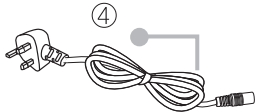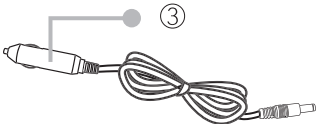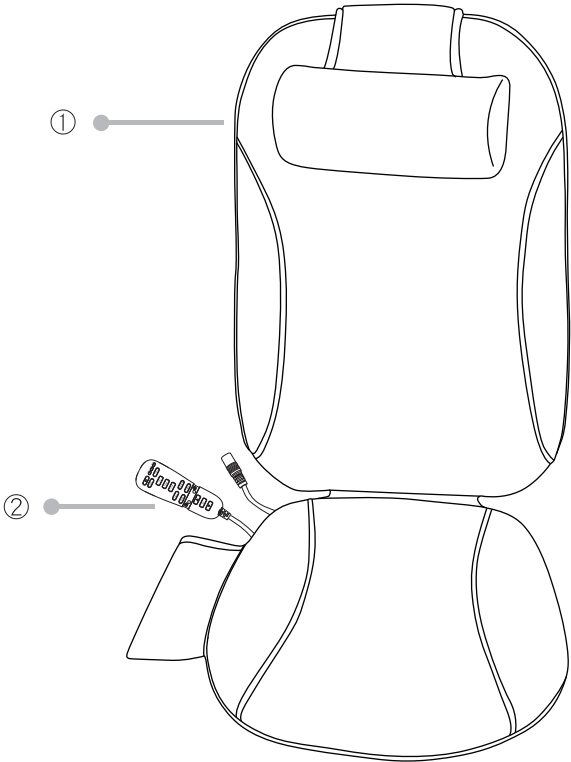Controller Diagram

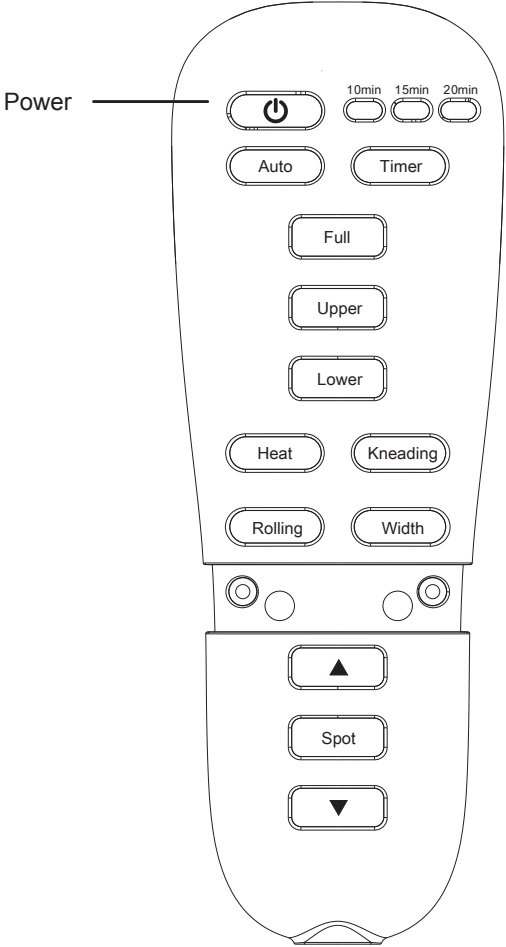Product Operation Guide

i. Placing the massager

- The massager equipped with a back strap that allows user to fasten onto most chairs.
- Place the massager on a chair, wrap the two ends of the straps round the back of the chair and secure them by the Velcro.

Operation instruction

i. [Power]

Insert the power plug into the electric outlet. Turn on the power. If the rollers are not at the bottom position, the light of the [Power] button will flash and the rollers will return to the bottom position.

Press the [Power] button, the product button [Auto] will flashing. The massager would automatically begin to run the full massage function with [Auto 1] mode.

Press again the [Auto] button, The product enter another [Auto 2] mode in which the [Full], [Heat], and [Rolling] button is activate and function.

The operation of massager for each usage is 15 minutes. After about 15 minutes of usage of the appliance, all the appliance function will stop.

Press the [Power] button, all the lights on the panel will turn off and all the buttons will stop functioning and the massager roller will set back to the original position.

Product Operation Guide

ii. [Rolling]

Press the [Rolling] button and the massager will perform the full back rolling, and heating process.

The [Full] back massage is a default for the [Rolling] mode but you can still select either [Upper] or [Lower] massage.

iii. [Kneading]

Press the [Kneading] button and the massager will perform the full back kneading and heating process.

The [Full] back massage is a default for the [Kneading] mode but you can still select either [Upper] or [Lower] massage.

iv. [Spot Up]/[Spot Down]

[Spot Up]/[Spot Down] function would work just when the [Kneading] function is operating.

The roller will fix to knead at a position when we press the [Spot Up]/[Spot Down] button.

The roller would move up to knead when we keep pressing the [Spot Up] button. It would fix until we stop pressing the [Spot Up] button.

Press the [Spot Down] button, and the roller would move down to knead and it would fix until we stop pressing the [Spot Down] button. Spotting function: the kneading function is only available when [Kneading] function is activated.

v. [Width]

After being activated, the [Width] function is only available when the automatic mode of [Kneading] is not enabled.

Press the [Width] button one time and the roller angle would be adjusted to another width.

Repeated press the [Width] button, the roller angle will adjust to another width.

Display

1. Kneading and Rolling massaging programs for full back, upper back or lower back relaxation.
2. Spot massage function for customized Kneading to relax the back area.
3. Heat function to increase the Kneading and Rolling effects and improve blood circulation.
4. Light and portable design.
5. Easy to use remote controller.
6. Auto shut-off after 15 minutes.
7. LED light indicators.

Effects by using the appliance

1. Relieves the muscle ache and pain.
2. Relieves the tension and reduces the stress.
3. Improves blood circulation.

When not in use, press the [Power] button, and please don't cut off the power immediately until the rollers have resumed to the bottom and the indicator light is off.

After use

Press the 'Power' button after use and the massager will stop operating. Massage rollers will resume to their original position. Unplug the appliance from the power outlet.

Product Maintenance

To store

- Keep the massager in the box or in a safe, dry and cool place.
- Keep the massager away from sharp edges or pointy objects that might cut or puncture the product's fabric surface.
- To avoid damage, do not wrap the power cord around the unit.
- Do not hang the unit with the controller cord.

To clean

- Avoid immersing unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Keep the appliance away from water or any liquid.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/ furniture polish or paint thinner to clean the product.

CAUTION

- Do not lean on the dual massage mechanism in the cushion when adjusting your body position to avoid tears on the products cushion surface.
- Do not jam or force any part of your body against the moving massage mechanism.
- Not recommended for use on wooden surfaces as the zipper may damage the wood.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.
- Driver is not recommended to use this product when driving as the heated function requires much more power. To avoid damaging your car's electricity system, connect the car plug to the product only after you activate your car's engine.
- Do not attempt to repair the massage cushion, please contact 'OGAWA Service Center' for any servicing required.

Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

Conditions	Possible Reason	Solution
Cant turn on	<ol style="list-style-type: none">1. Power not connect2. Remote not on	<ol style="list-style-type: none">1. Ensure the plug is properly connected to power source.2. Press the [Power] button
Sudden stop	<ol style="list-style-type: none">1. Power off2. Massager had over 15 minutes massage usage.3. Overheat protection4. Overload protection	<ol style="list-style-type: none">1. Check the power source2. Restart the massager3. Let the product cool, then restart the massager4. Restart the massager

Product Specifications

Product Name	Mobile Seat XE Duo
Model number of the Product	OZ0939
Power source	100-240V
Frequency	50-60Hz
Voltage	12V 2.5A
Power consumption	30W
Dimensions (Length/width/height)	L700xW179xH488 (mm)
Weight	5 KG

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F,YHC Tower,1Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel:+852 3465 6262
Fax:+852 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampine Industrial Avenue 5,
T5 @ Tampines
Singapore 528623

Tel:+065-6846 8588
Fax:+065-6846 7677

VIETNAM
OGAWA Vietnam Co Ltd
2nd floor,Petroland Tower,
12 Tan Trao,Tan Phu ward
District 7
Ho Chi Minh
Tel:+84 8 3775 3093/3775 3094
Fax:+84 8 3775 3095

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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HONG KONG, CHINA
 OGAWA Health Care International (HK) Ltd
 31/F,YHC Tower,1Sheung Yuet Road
 Kowloon Bay
 Kowloon
 Hong Kong
 Tel:+852 3465 6262
 Fax:+852 2410 0119

SINGAPORE
 OGAWA Health-Care Pte Ltd
 27 Tampine Industrial Avenue 5,
 T5 @ Tampines
 Singapore 528623

 Tel:+065-6846 8588
 Fax:+065-6846 7677

VIETNAM
 OGAWA Vietnam Co Ltd
 2nd floor,Petroland Tower,
 12 Tan Trao,Tan Phu ward
 District 7
 Ho Chi Minh
 Tel:+84 8 3775 3093/3775 3094
 Fax:+84 8 3775 3095

PHILIPPINES
 O-Healthcare Solution Phil.Inc.
 2303, 4th Floor,G&A Building,
 Chino Roces Avenue Extension,
 Makati City,Philippines
 Tel : +632 889 0676 / + 632 844 6851
 Fax: +632 812 0421

Ogawa is grateful
to have received
these awards:

Ogawa Wellness Hubs : Malaysia . Singapore . Vietnam . Philippines . Indonesia . Myanmar . India . China
Hong Kong . Taiwan . Turkey . Pakistan . Saudi Arabia . Mauritius . Russia . France
Canada