*mobile*seat **NX**

S-TRACK MASSAGE CUSHION

User Manual
OZ 0988 使用手册

*mobile*seat **NX**

S-TRACK MASSAGE CUSHION

User Manual

All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Controller Diagram	07
Product Operation Guide	08
Display	12
Product Maintenance	13
Troubleshooting	14
Product Specifications	14
Contact Us	15

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Do not use the appliance in the bathroom.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Place the appliance in a safe location.
- Do not use or sit on the appliance when driving.
- Do not expose the appliance to a heater or direct sunlight.

ii. Health precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease.
- Pregnancy.
- Any health concerns.

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Before sleeping.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

iv. Caution and warning

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by OGAWA.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to OGAWA Service Center for examination and repairs.
- Do not insert your hands or any foreign matters into the massager's vents and gaps.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- Do not lift or carry this appliance by the power cord or any cord.
- This product is only intended for personal use only. Commercial use is prohibited.

Safety Instructions

- Do not cover the appliance during operation. As this may result in malfunction or fire.
- Do not crush the appliance and avoid any sharp folds.

V. Note


- This appliance can be used by children at the age of 8 or above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under guidance or if they have been informed of the risks involved and instructions concerning use of the appliance in a safe way.
- Children shall not play on the appliance.
- Cleaning and maintenance shall not be performed by children without supervision.
- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- Extensive usage could lead to excessive heating and a shorter product lifespan. Should this occur, discontinue the use and allow the unit to cool before operating again.
- It is normal to hear the noises during use, because they are resulted from vibration of product.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.
- If damaged, the supply cord must be repaired by the manufacturer, after-sales service provider or qualified persons in similar field so as to avoid a hazard.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Warranty Claims


1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

- A) Back Cushion
- B) Detachable Backrest Flap
- C) Seat Cushion
- D) Cable plug
- E) Controller

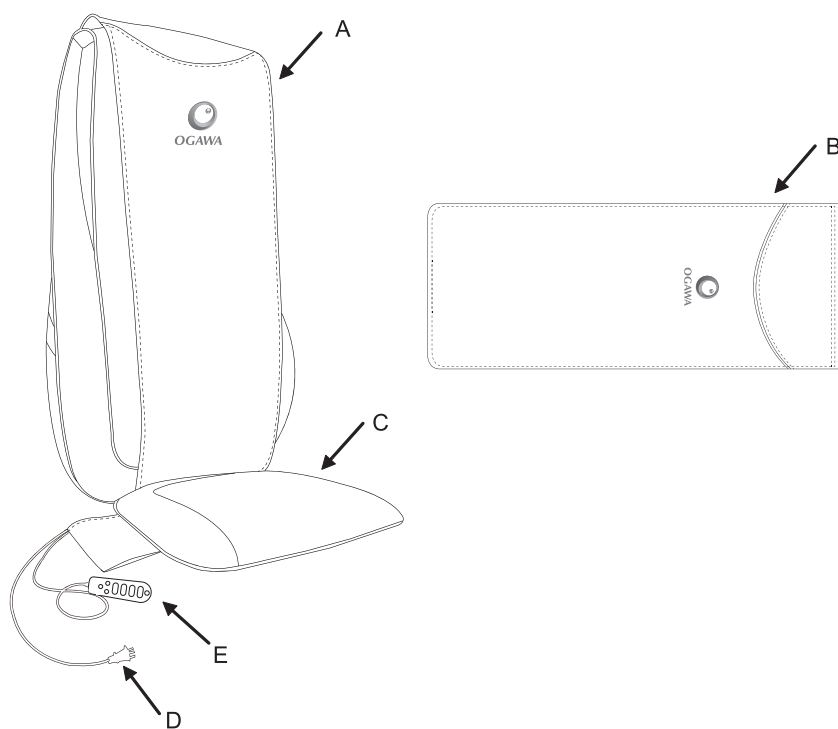Controller Diagram

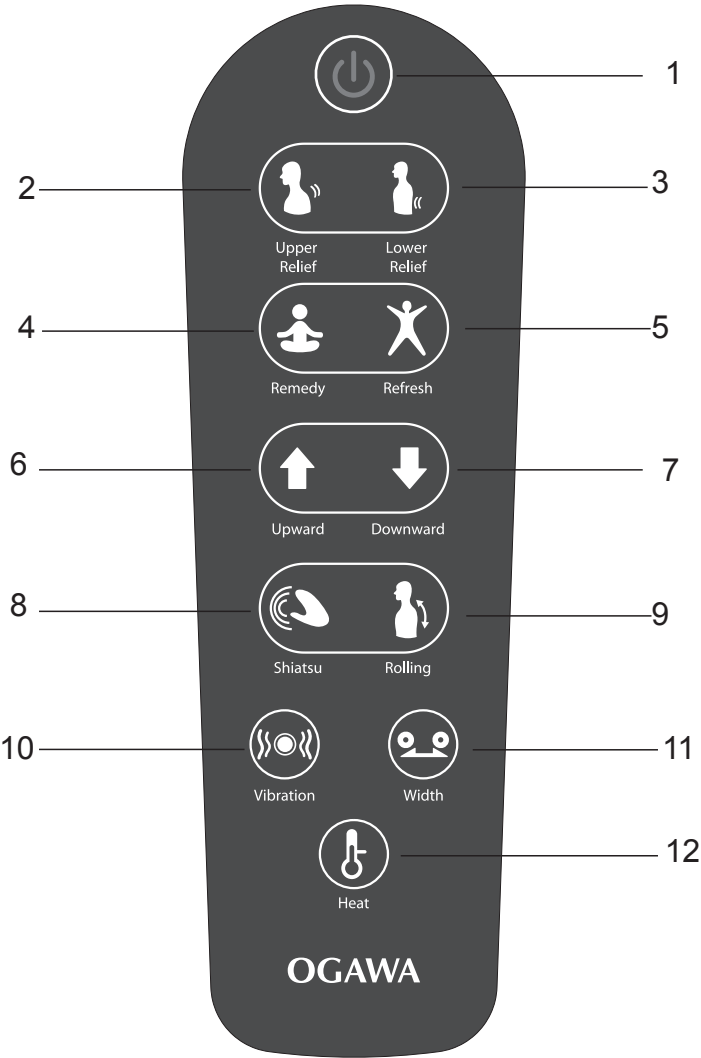Product Operation Guide

i. Preparation before use

- A screw is applied on the back of the product to protect the massager apparatus during shipment. Remove the screw before use.
- Use the equipped wrench to remove the screw and insert the power cord provided into the opening.

ii. Placing the massager

- The massager is equipped with back strap that allow user to faster onto most chairs.
- Place the massager on a chair, wrap the two end of the strap at the back of the chair then secure with Velcro.

iii. Backrest flap

- A detachable backrest flap is provided for a gentler massage.
- Remove the backrest flap if you prefer more intense massage.

Product Operation Guide

1. Power
Insert the power plug into the electric outlet. Turn on the power.
If the roller are not at the bottom position, the light of the Power button will
2. Upper Relief
This key is valid when Power is on.
Press this button to activate the Upper Relief function, the indicator light turn on.
Press this button again to close the Upper Relief function, the indicator light turn off.
3. Lower Relief
This key is valid when Power is on.
Press this button to activate the Lower Relief function, the indicator light turn on.
Press this button again to close the Lower Relief Function, the indicator light turn off.
4. Remedy
This key is valid when Power is on.
Press this button to activate the Remedy function, the indicator light turn on.
Press this button again to close the Remedy function, the indicator light turn off.
5. Refresh
This key is valid when Power is on.
Press this button to activate the Refresh function, the indicator light turn on.
Press this button again to close the Refresh function, the indicator light turn off.
6. Spot upward
This key is valid in the spot / upper relief / lower relief / remedy / refresh. Long press this button, the roller move upward, the indicator light flashes; if the roller touch the limit sensor, roller will stop, the indicator light turn on; release the button, the roller will stop moving, indicator light turn off. At this point roller will change to kneading.
7. Spot downward
This key is valid in the spot / upper relief / lower relief / remedy / refresh. Long press this button, the roller move upward, the indicator light flashes; if the roller touch the limit sensor, roller will stop, the indicator light turn on; release the button, the roller will stop moving, indicator light turn off. At this point roller will change to kneading.
8. Shiatsu
This key is valid when Power is on.
Press this button, turn on (shiatsu + rolling) slow mode., the indicator light turn on.
Press this button second time, turn on (shiatsu + rolling) fast mode, the indicator light turn on.
Press this button again, turn off (shiatsu + rolling) function, the indicator light turn off.

Product Operation Guide

9. Rolling

This key is valid when Power is on.

Press this button, turn on rolling in slow mode.,the indicator light turn on.

Press this button second time, turn on rolling in fast mode, the indicator ligh turn on.

Press this button again, turn off rolling) function, the indicator light turn off.

10. Vibration

This key is valid when Power is on.

Press this button to activate the vibration, the indicator light turn on.

Press this button again to close the vibration, the indicator light turn off

11. Width adjustment

This key is valid when rolling function is on.

Press this button to adjust the massage roller width.

12. Heat

This key is valid in the spot / Upper Relief / Lower Relief / Remedy / Refresh

Press this button to activate the heating, the indicator light turn on.

Press this button again to close the heating, the indicator light turn off

Display

1. Shiatsu and rolling massaging programs can relax your back.
2. Fixed-position massaging function provides Shiatsu of specific position.
3. Seat cushion vibration relieves the lower body tension.
4. Heat function can increase the Shiatsu and Rolling effect and improve blood circulation.
5. Light and portable design.
6. Easy to use remote controller.
7. Auto shut-off after 15 minutes.
8. LED light indicators.

Effects by using the appliance

1. Relieves the muscle ache and pain.
2. Relieves the tension and reduce the stress.
3. Improves blood circulation.

When you do not use the product, press the [Power] button, and please don't cut off the power immediately until the rollers have resumed to the bottom and the indicator light is off.

Press the [Power] button after use and the massager will stop operating. Massage roller will resume to their original position. Unplug the appliance from power outlet.

Product Maintenance

TO STORE

- Keep massager in the box or a safe, dry and cool place.
- Keep away the massager from sharp edges or pointy objects that might cut or puncture the product's fabric surface.
- Do not wrap the power cord around the unit.
- Do not hang the unit with the controller cord.

TO CLEAN

- Avoid immersing the unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Keep the appliance away from water or any liquid.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish or paint thinner to clean the product.

CAUTION

- Do not lean on the dual massage mechanism in the cushion when adjusting your body position to avoid tears on the products cushion surface.
- Do not jam or force any part of your body against the moving massage mechanism.
- Not recommended for use on wooden surfaces as the zipper may damage the wood.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.
- Do not attempt to repair the product, please contact 'OGAWA Service Center' for any servicing required.
- Beware when put on light color PU / fabric seat, some material may have color dyeing
- Beware when put printing document on product, some material may have color dyeing

Troubleshooting

Symptoms	Solutions
Not operating.	<ul style="list-style-type: none"> • Use the product when electricity is restored. • Reconnect the power cord. • Turn on the power switch.
No light & heating after pressing the heat switch.	<ul style="list-style-type: none"> • Press the 'Shiatsu' switch before turning the 'Heat' switch. • Press the 'Rolling' switch before turning on the 'Heat' switch.
During operation, the controller is working but the massager is not working.	<ul style="list-style-type: none"> • The massager could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and stop using the massager for 30 to 50 minutes to facilitate heat dissipation of the engine.
Back massage function does not work.	<ul style="list-style-type: none"> • Turn on the power switch. • Reconnect the power cord.
Seat vibration does not work.	<ul style="list-style-type: none"> • Turn on the 'Vibration' button. • Reconnect the power cord.
Sudden stop.	<ul style="list-style-type: none"> • Use the product when electricity is restored. • Reconnect the power cord.

Product Specifications

Product Name	Mobile Seat NX
Model number of the product	OZ 0988
Power source	220-240V~
Frequency	50~60Hz
Output Voltage	24V --- 2.5A
Power consumption	60W
Dimension (Length/width/height)	L810 x W480 x H210 (mm)
Weight	7.87kg

Contact Us

MALAYSIA
Healthy World Lifestyle Sdn Bhd
No. 22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel : +603-5121 4286
Fax: +603-5121 4386

HONG KONG, CHINA
OGAWA Health Care International (HK) Ltd
31/F, YHC Tower, 1 Sheung Yuet Road
Kowloon Bay
Kowloon
Hong Kong
Tel :+852 3465 6262
Fax:+852 2410 0119

SINGAPORE
OGAWA Health-Care Pte Ltd
27 Tampines Industrial Avenue 5,
T5 @ Tampines
Singapore 528623
Tel : +065-6846 8588
Fax: +065-6846 7677

VIETNAM
OGAWA Vietnam Co Ltd
2nd floor, Petroland Tower,
12 Tan Trao, Tan Phu ward
District 7
Ho Chi Minh
Tel : +84 8 3775 3093 / 3775 3094
Fax: +84 8 3775 3095

INDONESIA
PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel : +62 21-6983 1123
Fax: +62 21-6983 1124

PHILIPPINES
O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel : +632 889 0676 / +632 844 6851
Fax: +632 812 0421

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Thailand . Korea . Australia . Dubai . K.S.A . Nigeria . Egypt . Lebanon . Denmark
Sweden . U.S.A